

rose plastic medical packaging

General catalog





Safe packaging solutions

for medical devices, dental, healthcare and laboratory

Packaging in the field of medical devices and healthcare today has to meet particularly strict requirements. The spectrum ranges from protection against mechanical damage and effective protection against contamination to double-sterile-barrier packaging for maximum protection against the ingress of germs.

rose plastic medical packaging specializes precisely in this. We not only develop and produce all our hard plastic packaging products ourselves, but also all the necessary molds.

We also develop tailor-made packaging solutions in close consultation with our customers.

So whenever you need high-quality, practical and secure protective packaging, you've come to the right place.

We will be happy to advise you on the many possible uses and advantages of our products.

Best wishes

Bernd Gerber

Managing Director
rose plastic medical packaging GmbH

2013 - 2023
10 years



We are your flexible partner

International. Innovative. Ideal.

rose plastic medical packaging specializes in the development and manufacturing of packaging for medical devices such as instruments, implants and accessories, as well as for healthcare, laboratory and diagnostic products.

With production facilities in Germany, the USA, Brazil, China and India, as well as sales offices worldwide, the company offers cleanroom quality in accordance with ISO 8 and ISO 9 as well as reliable service.


rose plastic medical packaging is a company of the rose plastic group. rose plastic is the world market leader in plastic packaging for precision tools. The company, headquartered in Southern Germany, has been known since 1953 for the development and production of innovative plastic packaging and its outstanding know-how in injection molding and blow molding.

All from one source

From idea to finished package



Because we have all the essential process steps in our own hands, we have full control. This high vertical range of manufacture guarantees our customers a consistently high product and process quality.

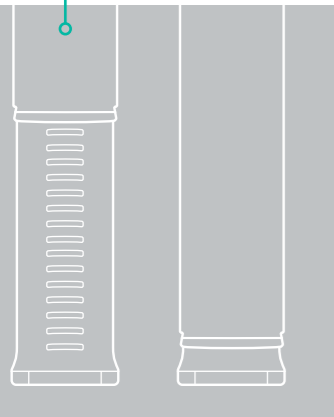

1 The idea

DESIGN AND DEVELOPMENT

A good product is first and foremost based on a good idea that not only takes the customer's wishes into account, but also often exceeds them. The experts in our highly efficient in-house development department take the design and construction of our packaging into their own hands right from the start. They and the proximity to our customers are the reasons for our well-designed, high-quality products.



2 The tool

TOOL SHOP

We produce our packaging using high-precision blow molding and injection molding tools. But who produces our tools? We also do this ourselves – from the design right through to CNC-supported production. Because decades of experience have taught us the exact requirements for our production tools. Through the years, we have developed and integrated numerous special features that ensure superior process and product quality.

3 The tests

TECHNICAL CENTER

At our in-house technical center we put new tools to the acid test, carry out exhaustive production tests, and give valuable feedback to our development team and tool shop. Our colleagues in this department also perform exacting material tests and trials.



4 The series

PRODUCTION

Our packaging products only go into series production (injection molding and blow molding) once they have been fine-tuned and have successfully passed various tests and inspections. We have also introduced numerous special features into our production process. These allow us to work with higher process reliability in order to achieve higher product quality, while ensuring economical production.


5 The logistics

STORAGE AND SHIPPING

Our warehouses, which are designed for maximum capacity, and our experienced employees ensure that our products are well packed and reach you as soon as possible. Which in turn enables you to quickly deliver well-packaged products to your own customers.



Where is our packaging in use?

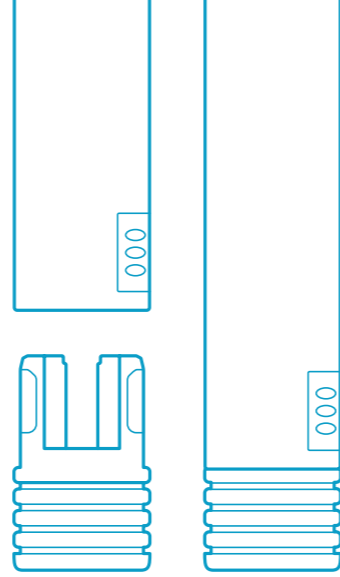

Medical device

GREAT THINGS DESERVE EXCELLENT PACKAGING

Whether you need packaging for a delicate implant or surgical and medical instruments - rose plastic medical packaging can supply both primary and secondary packaging for a huge range of medical devices.



Dental

NO DRILL IS TOO SMALL FOR OUR PACKAGING

Whether you need packaging for dentistry supplies or dental technology, dental surgery or implantology products, all our packaging has one thing in common - it always meets the highest safety and quality standards. And the possible applications are incredibly diverse: Dental instruments, inlays, crowns, handpieces and contra-angle handpieces, drills and files - our well-thought-out packaging provides a safe home for almost everything.



Healthcare

THE PERFECT PACKAGING FOR YOUR HEALTHCARE PRODUCT

Even the best healthcare product stands or falls by the right packaging.. Here you will find a variety of packaging solutions that not only protect your product perfectly, but also present it in an appealing way.



Laboratories and diagnostics

NEATLY STOWED AND ALWAYS READY TO HAND

We help you keep order in the lab and workplace. We supply intelligent packaging solutions for the lab and diagnosis. Our range also includes the AQRADO® modular storage system.





Precise. Clean. Reliable.

In order to deliver the quality you need, we manufacture our products in state-of-the-art cleanrooms. We also ensure consistently high standards via our certified quality management system in accordance with EN ISO 13485.

Protective and transport packaging for sterile and non-sterile applications

To ensure we can supply clean packaging systems for your products, we manufacture within a strictly controlled environment that complies with cleanroom classification ISO Class 9 (EN ISO 14644-1) and ISO Class 8. In addition, our comprehensive standard product range is manufactured within the framework of EN ISO 9001.

Our cleanrooms are located at our headquarters in Hergensweiler (Germany) and in California (Pennsylvania, USA). We also offer production in a clean environment at our location in Sorocaba (Brazil). Naturally, we use medical grade materials to meet our customers' exacting demands.

These include, in particular, common materials such as PETG, PE, PP and TPE.

We deliver certified quality – guaranteed!

At our location in Hergensweiler (Germany), rose plastic medical packaging is certified in accordance with the following international standards:



innovations in plastic packaging

PRODUCT RANGE

Protective Cover



TipCover	14
UniCover	16

Tubes



BlockPack	18
BlockPack TabletPunch	20
DentalImplantPack	22
TwistPack	24
FixInTube®	26
ScrewPack	28
MicroPack XPress	30
PedicleScrew Pack	32
PlugPack	34
QuadroPack	36
QuadroPack TipLine	38
RectangularPack	40
RectangularPack TipLine	43
TelePack	44
TipTube	46
TopPack XPress	48
TubeInTube®	52
TwistPack Plus	54

Boxes



ConsumerBox	56
SplitPack	58
TapePack	61
UniBox	62

Cases



RoseCase ErgoLine	64
RoseCase ProTec	66

Accessory




AQURADO®	68
----------------	----





Tip Cover

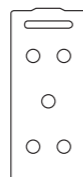

Optimum protection for delicate tips and edges

- ideal for the protection of tips and edges of surgical instruments
- for sterile and non-sterile applications
- secure loading and removal
- enhanced tip-protection
- security and protection for your product
- easy handling
- safe transportation
- flexibly adjusts to the product geometry
- 6 sizes available for various diameter ranges



The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.

Product Range

Article	Outer dimensions W x H x D	Suitable for	PU
TCV0 0.5x28/1	10x28x5mm	0,5 - 2,5mm	1000
TCV1 2.5x28/1	9,3x28,3x4,8mm	2,5 - 4,0mm	1000
TCV2 4x28/1	12,5x28,4x7,5mm	4,0 - 7,0mm	1000
TCV3 7x28/1	17,0x28,5x12mm	7,0 - 10,5mm	1000
TCV4 10.5x28/1	23,3x28,7x16mm	10,5 - 14,5mm	500
TCV5 14.5x36/1	29,6x35,8x22mm	14,5 - 18,0mm	250



MATERIAL

- TPE

PRODUCTION ENVIRONMENT

- cleanroom production

STANDARD COLOR



“It is the small things in life which count.”
Charles Comiskey



UniCover




Flexible protective cover for different applications

- abrasion and puncture resistant
- great protection of the device surface
- suitable to gamma, ethylene oxide and eBeam sterilization
- protects sharp and sensitive medical instruments, implants, screws etc.
- UniCover is made of TPU (thermoplastic polyurethane). This material is characterized by excellent abrasion resistance and is suitable for various sterilization processes.
- intuitive handling thanks to "one-handed removal" feature also helps healthcare professionals save time when using it
- When closed, the packaged product is held inside the packaging by the closure studs. By applying lateral pressure in the area around these studs, the packaging is compressed, allowing the product to be removed.
- available with a cap on request
- internal geometry can be customized




MATERIAL
- Thermoplastic polyurethane (TPU)

LENGTH

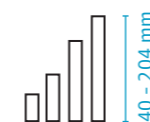

WIDTH

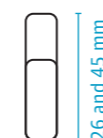

PRODUCTION ENVIRONMENT

- cleanroom production

STANDARD COLOR



Product Range

Article	Width (mm)	Length (mm)	Depth (mm)	Internal geometry
UCV 26x40x8/1	26	40	8	4 closure studs
UCV 26x40x8/2	26	40	8	-
UCV 26x98x8/1	26	98	8	4 closure studs
UCV 26x98x8/2	26	98	8	-
UCV 45x134x10/1	45	134	11	4 closure studs
UCV 45x134x10/2	45	134	11	-
UCV 45x204x10/1	45	204	11	4 closure studs
UCV 45x204x10/2	45	204	11	-


All dimensions are external dimensions.



UniCover sets new standards in terms of safety and user-friendliness.



BlockPack

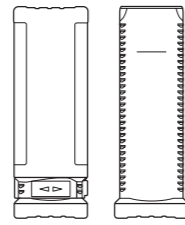

Square telescopic packaging tubes with twist lock closure feature and ratchet length adjustment

- highly secure bayonet-type twist lock closure, requires minimal force in opening and closing the pack
- contoured tube ends and lateral support grooves are integrated together to provide a cushioning effect
- with ratchet length adjustment
- ergonomic product design offer shape stability
- speedy handling
- safe transportation
- excellent for stacking and storage
- durable and reusable
- can be used as primary or secondary packaging



The ultimate, universal packaging tube with telescopic length adjustment – for maximum product protection

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.



PACKAGING
Outer parts and inner parts are packed separately




MATERIAL
— Polyethylene (PE)




COLOR
— PE natural

USABLE LENGTHS

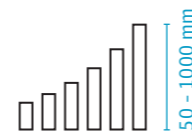

Ø DIAMETER

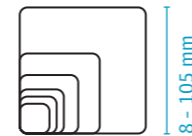

PRODUCTION ENVIRONMENT

— various sizes available in cleanroom production

STANDARD COLOR



Product Range

Article	Ø	Usable length	PU	Article	Ø	Usable length	PU
BK 8 050	8mm	50 - 80mm	3000	BK 32 060	32mm	60 - 80mm	800
BK 8 080	8mm	80 - 120mm	3000	BK 32 080	32mm	80 - 120mm	600
BK 8 120	8mm	120 - 200mm	1500	BK 32 120	32mm	120 - 200mm	400
BK 8 200	8mm	200 - 350mm	2000	BK 32 200	32mm	200 - 350mm	360
BK 10 050	10mm	50 - 80mm	2100	BK 32 350	32mm	350 - 620mm	140
BK 10 080	10mm	80 - 120mm	2000	BK 35 050	35mm	50 - 70mm	900
BK 10 120	10mm	120 - 200mm	1170	BK 35 060	35mm	60 - 80mm	700
BK 10 200	10mm	200 - 350mm	1100	BK 35 080	35mm	80 - 120mm	500
BK 12 050	12mm	50 - 80mm	2200	BK 35 120	35mm	120 - 200mm	300
BK 12 080	12mm	80 - 120mm	1550	BK 35 200	35mm	200 - 350mm	300
BK 12 120	12mm	120 - 200mm	800	BK 35 350	35mm	350 - 620mm	159
BK 12 200	12mm	200 - 350mm	800	BK 40 080	40mm	80 - 120mm	450
BK 12 250	12mm	250 - 460mm	700	BK 40 120	40mm	120 - 200mm	300
BK 14 050	14mm	50 - 80mm	1700	BK 40 200	40mm	200 - 350mm	240
BK 14 080	14mm	80 - 120mm	1100	BK 40 350	40mm	350 - 620mm	120
BK 14 120	14mm	120 - 200mm	800	BK 45 060	45mm	60 - 80mm	450
BK 14 200	14mm	200 - 350mm	600	BK 45 080	45mm	80 - 120mm	350
BK 16 050	16mm	50 - 80mm	1600	BK 45 120	45mm	120 - 200mm	275
BK 16 080	16mm	80 - 120mm	850	BK 45 200	45mm	200 - 350mm	160
BK 16 120	16mm	120 - 200mm	1200	BK 45 350	45mm	350 - 620mm	90
BK 16 200	16mm	200 - 350mm	1200	BK 55 060	55mm	60 - 80mm	250
BK 18 050	18mm	50 - 80mm	1300	BK 55 080	55mm	80 - 120mm	200
BK 18 080	18mm	80 - 120mm	850	BK 55 120	55mm	120 - 200mm	200
BK 18 120	18mm	120 - 200mm	1050	BK 55 200	55mm	200 - 350mm	100
BK 18 200	18mm	200 - 350mm	400	BK 55 350	55mm	350 - 620mm	70
BK 20 050	20mm	50 - 80mm	1200	BK 65 080	65mm	80 - 120mm	210
BK 20 080	20mm	80 - 120mm	1300	BK 65 120	65mm	120 - 200mm	120
BK 20 120	20mm	120 - 200mm	900	BK 65 200	65mm	200 - 350mm	90
BK 20 200	20mm	200 - 350mm	800	BK 65 350	65mm	350 - 620mm	48
BK 20 350	20mm	350 - 620mm	500	BK 75 080	75mm	80 - 120mm	168
BK 22 050	22mm	50 - 80mm	1000	BK 75 120	75mm	120 - 200mm	96
BK 22 080	22mm	80 - 120mm	1200	BK 75 200	75mm	200 - 350mm	58
BK 22 120	22mm	120 - 200mm	800	BK 85 080	85mm	80 - 120mm	90
BK 22 200	22mm	200 - 350mm	600	BK 85 120	85mm	120 - 200mm	72
BK 24 050	24mm	50 - 80mm	750	BK 85 200	85mm	200 - 350mm	48
BK 24 080	24mm	80 - 120mm	1000	BK 85 350	85mm	350 - 620mm	30
BK 24 120	24mm	120 - 200mm	700	BK 85 550	85mm	550 - 1000mm	20
BK 24 200	24mm	200 - 350mm	550	BK 95 080	95mm	80 - 120mm	80
BK 26 050	26mm	50 - 60mm	700	BK 95 120	95mm	120 - 200mm	60
BK 26 060	26mm	60 - 80mm	700	BK 95 200	95mm	200 - 350mm	30
BK 26 080	26mm	80 - 120mm	850	BK 95 350	95mm	350 - 620mm	20
BK 26 120	26mm	120 - 200mm	600	BK 105 080	105mm	80 - 120mm	80
BK 26 200	26mm	200 - 350mm	450	BK 105 120	105mm	120 - 200mm	48
BK 26 350	26mm	350 - 620mm	288	BK 105 200	105mm	200 - 350mm	30
BK 32 050	32mm	50 - 70mm	1200	BK 105 350	105mm	350 - 620mm	20

Awards



BlockPack TabletPunch



Packaging solution for tableting punches with twist lock closure feature

- the sensitive tip of the tableting punch is suspended for optimum protection
- highly secure bayonet-type twist lock closure, requires minimal force in opening and closing the pack
- contoured tube ends and lateral grooves provide protection for the punch
- ergonomic product design offers shape stability
- easy handling
- excellent for stacking and storage
- durable and re-usable
- two sizes accommodate the most common tableting tool designs
- please note that rose plastic BlockPack TabletPunch must not be sold or resold in the UK



The perfect packaging solution for optimal protection of high-quality tableting punches



PACKAGING
Outer parts and inner parts are packed separately



MATERIAL
- Polyvinylchloride (PVC)




COLOR
- PVC clear

USABLE LENGTHS



Ø DIAMETER

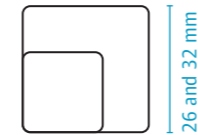

PRODUCTION ENVIRONMENT

- Industrial standard



STANDARD COLOR




Product Range

Article	Inner diam.	Shank diam.	Usable length fix	PU
BKT 26 140 A+I	26mm	19.5mm	140mm	600
BKT 32 140 A+I	32mm	26mm	140mm	400

A = outer part; I = inner part

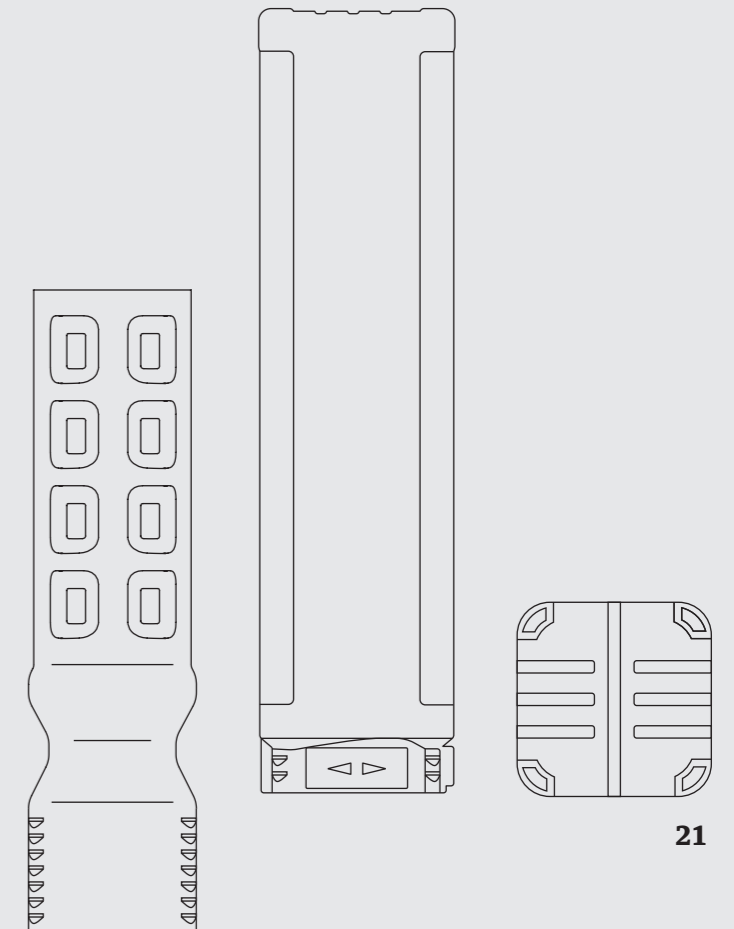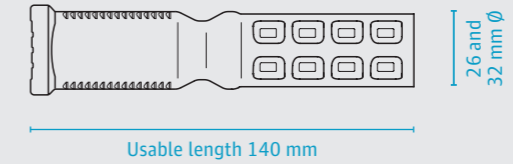

Dental ImplantPack

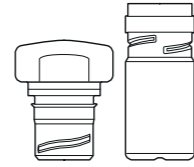

Primary packaging for dental implants

- high quality single packaging for the protection of dental implants with and without implant mount
- the implant is clearly visible in the highly transparent packaging
- the implant holder fixes the implant and the cover screw firmly in the packaging
- simple handling adapted to the operation process
- easy and quick removal of the implant during the operation by the dentist
- by turning the cover on the implant holder, the implant lays free and can be directly picked up for insertion
- suitable for dental implants with \varnothing 2.8 - 6.5 mm and lengths up to 22 mm
- matching titanium inserts or the clamping on the implant mount prevent contact between plastic and implant
- the different diameters of the implants can be covered via a titanium ring, the titanium length stop can be flexibly adjusted to the implant length
- the rollaway protection integrated in the closure provides additional safety in the sterile operating area
- sterile barrier system possible with additional blister with Tyvek® or sterilization pouch



Easy handling, safe removal and flexible implant holding: the Dental ImplantPack fully convinces.

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.



PACKAGING
Outer parts and inner parts are packed separately



MATERIAL


- Tube: PETG
- Closure: Polyethylene (PE)/TPE
- Implant holder: Polyethylene (PE)
- Cover: PETG



COLOR

- Tube: PETG clear
- Closure: PE/TPE natural

USABLE LENGTHS

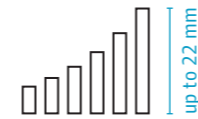

Ø DIAMETER



PRODUCTION ENVIRONMENT

- Cleanroom production


STANDARD COLOR



Product Range

Article	Implant \varnothing max.	Implant length max.	Suitable for implants	Cover screw holding	Implant holding
DIP 16/1	6.5mm	22mm	without implant mount	universal (adaptable)	Titanium ring and plate / titanium tube (not included)
DIP 16/2	6.5mm	18mm	with implant mount	universal (adaptable)	Clamping at implant mount (adaptable)

TwistPack

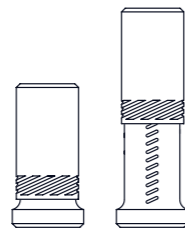

Round telescopic protective packaging tubes with twist lock mechanism and roll-away protection

- twist lock mechanism holds contents securely in place
- infinitely variable telescopic length adjustment
- hexagon tube bottom prevents rolling
- security and protection for your product
- speedy efficient packaging
- speedy handling
- safe transportation
- durable and reusable
- can be used as primary or secondary packaging



PACKAGING

Outer parts and inner parts are packed separately; DP 65 350, DP 80 and DP 105 outer parts and inner parts are packed assembled



MATERIAL


- Polypropylen (PP)



COLOR

- PP transparent

USABLE LENGTHS

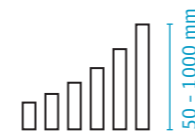

Ø DIAMETER



PRODUCTION ENVIRONMENT

- various sizes available in cleanroom production

STANDARD COLOR



Product Range

Article	Inner diam.	Usable length	PU	Article	Inner diam.	Usable length	PU
DP 10 050	10mm	50 - 70mm	2100	DP 35 080	35mm	80 - 120mm	600
DP 10 080	10mm	80 - 120mm	2000	DP 35 120	35mm	120 - 200mm	400
DP 10 120	10mm	120 - 200mm	2100	DP 35 200	35mm	200 - 350mm	310
DP 13 050	13mm	50 - 70mm	2200	DP 35 350	35mm	350 - 620mm	180
DP 13 080	13mm	80 - 120mm	1300	DP 40 050	40mm	50 - 60mm	350
DP 13 120	13mm	120 - 200mm	1600	DP 40 060	40mm	60 - 80mm	350
DP 13 200	13mm	200 - 350mm	1500	DP 40 080	40mm	80 - 120mm	500
DP 16 050	16mm	50 - 70mm	1700	DP 40 120	40mm	120 - 200mm	300
DP 16 080	16mm	80 - 120mm	1000	DP 40 200	40mm	200 - 350mm	250
DP 16 120	16mm	120 - 200mm	1200	DP 40 250	40mm	250 - 400mm	210
DP 16 200	16mm	200 - 350mm	1000	DP 40 350	40mm	350 - 620mm	150
DP 16 350	16mm	350 - 620mm	600	DP 45 050	45mm	50 - 60mm	350
DP 16 410	16mm	410 - 770mm	500	DP 45 060	45mm	60 - 80mm	600
DP 20 050	20mm	50 - 70mm	1500	DP 45 080	45mm	80 - 120mm	400
DP 20 080	20mm	80 - 120mm	1700	DP 45 120	45mm	120 - 200mm	250
DP 20 120	20mm	120 - 200mm	1000	DP 45 200	45mm	200 - 350mm	200
DP 20 200	20mm	200 - 350mm	900	DP 45 350	45mm	350 - 620mm	115
DP 22 050	22mm	50 - 60mm	1100	DP 50 050	50mm	50 - 60mm	600
DP 22 080	22mm	80 - 120mm	1400	DP 50 060	50mm	60 - 80mm	500
DP 22 120	22mm	120 - 200mm	850	DP 50 080	50mm	80 - 120mm	330
DP 22 200	22mm	200 - 350mm	700	DP 50 120	50mm	120 - 200mm	220
DP 22 350	22mm	350 - 620mm	390	DP 50 200	50mm	200 - 350mm	150
DP 22 410	22mm	410 - 770mm	300	DP 50 350	50mm	350 - 620mm	95
DP 24 050	24mm	50 - 60mm	900	DP 55 050	55mm	50 - 60mm	500
DP 24 080	24mm	80 - 120mm	1100	DP 55 060	55mm	60 - 80mm	400
DP 24 120	24mm	120 - 200mm	750	DP 55 080	55mm	80 - 120mm	300
DP 24 200	24mm	200 - 350mm	600	DP 55 120	55mm	120 - 200mm	190
DP 26 050	26mm	50 - 60mm	800	DP 55 200	55mm	200 - 350mm	135
DP 26 080	26mm	80 - 120mm	1000	DP 55 350	55mm	350 - 620mm	70
DP 26 120	26mm	120 - 200mm	600	DP 65 050	65mm	50 - 60mm	350
DP 26 200	26mm	200 - 350mm	450	DP 65 060	65mm	60 - 80mm	300
DP 26 250	26mm	250 - 400mm	420	DP 65 080	65mm	80 - 120mm	220
DP 26 350	26mm	350 - 620mm	280	DP 65 120	65mm	120 - 200mm	150
DP 26 410	26mm	410 - 770mm	240	DP 65 200	65mm	200 - 350mm	96
DP 28 050	28mm	50 - 60mm	800	DP 65 350	65mm	350 - 620mm	46
DP 28 080	28mm	80 - 120mm	950	DP 80 060	85mm	60 - 80mm	100
DP 28 120	28mm	120 - 200mm	600	DP 80 080	85mm	80 - 100mm	108
DP 28 200	28mm	200 - 350mm	450	DP 80 120	85mm	120 - 180mm	72
DP 30 050	30mm	50 - 60mm	700	DP 80 200	85mm	200 - 350mm	50
DP 30 080	30mm	80 - 120mm	800	DP 80 350	85mm	350 - 550mm	18
DP 30 120	30mm	120 - 200mm	500	DP 80 550	85mm	550 - 1000mm	18
DP 30 200	30mm	200 - 350mm	400	DP 105 080	105mm	80 - 100mm	72
DP 30 350	30mm	350 - 620mm	230	DP 105 120	105mm	120 - 180mm	52
DP 30 410	30mm	410 - 770mm	200	DP 105 200	105mm	200 - 350mm	30
DP 35 050	35mm	50 - 60mm	500				
DP 35 060	35mm	60 - 80mm	900				



As one of the world's best-selling packaging tubes, the TwistPack excels thanks to practical features and almost limitless applications

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.

FixInTube®

Primary packaging for headed implant screws

- high quality individual packaging for implants offering excellent protection
- the special inner shape of the tube and the TPE/PE closure fix the implant firmly in the packaging
- contactless and safe removal of the implant screw using an attaching tool or surgical driver
- optimum security and protection for product and users
- high clarity packaging for perfect product visibility
- sterile barrier system possible by additional blister with Tyvek® or Pouch



FixInTube® provides a clear view. In the highly transparent tube, the contents are crystal clear visible at all times.

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.



PACKAGING
Outer parts and inner parts are packed separately




MATERIAL
— Tube: PETG
— Closure: combination of Polyethylene (PE) and TPE



COLOR
— Tube: PETG clear
— Closure: PE/TPE natural

USABLE LENGTHS

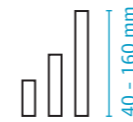

Ø DIAMETER



PRODUCTION ENVIRONMENT

— Cleanroom production


STANDARD COLOR



Product Range

Article	Screw thread Ø	Screw length max.	Diameter	Overall length
FT 10 040-2.0/1	2mm	40mm	10mm	77mm
FT 10 040-2.5/1	2.5mm	40mm	10mm	77mm
FT 10 040-3.0/1	3mm	40mm	10mm	77mm
FT 10 040-3.5/1	3.5mm	40mm	10mm	77mm
FT 10 040-6.0/1	6mm	40mm	10mm	77mm
FT 10 080-4.0/1	4mm	80mm	10mm	177mm
FT 10 080-5.0/1	5mm	80mm	10mm	177mm
FT 10 080-6.0/1	6mm	80mm	10mm	177mm
FT 10 160-7.0/1	7mm	160mm	10mm	197mm
FT 10 160-8.0/1	8mm	160mm	10mm	197mm

ScrewPack

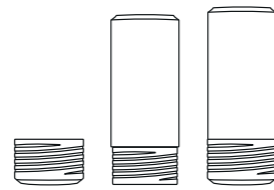

Round protective packaging tubes with secure screw closure

- length is produced to customer requirements
- universal use for products of different dimensions
- secure screw closure
- thread is positioned to customer requirements
- security and protection for your product
- speedy efficient packaging
- speedy handling
- safe transportation
- durable and reusable
- can be used as primary or secondary packaging



PACKAGING
Outer parts and inner parts are packed separately




MATERIAL
— Polypropylene (PP)



COLOR
— PP transparent

USABLE LENGTHS

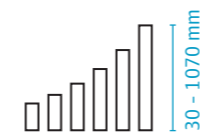

Ø DIAMETER



PRODUCTION ENVIRONMENT


- various sizes available in cleanroom production

STANDARD COLOR



Product Range


Article	Inner diam.	Length outer part min.- max.	Length inner part min.- max.	Total usable length
GP 25	28mm	25-407mm	5-387mm	794mm
GP 35	37mm	25-407mm	5-387mm	794mm
GP 40	41mm	35-390mm	5-360mm	750mm
GP 45	46mm	35-715mm	5-685mm	1400mm
GP 55	56mm	35-575mm	5-545mm	1020mm
GP 65	66mm	45-520mm	5-480mm	1000mm
GP 75	76mm	55-525mm	5-480mm	1005mm
GP 80	79mm	55-540mm	5-490mm	1030mm
GP 85	85mm	55mm	35mm	90mm
GP 90	91mm	55-485mm	5-436mm	921mm
GP 100	98mm	55-540mm	5-490mm	1030mm
GP 105	107mm	55-540mm	5-490mm	1030mm
GP 125	125mm	70-585mm	10-525mm	1110mm
GP 160	160mm	110-575mm	30-495mm	1070mm

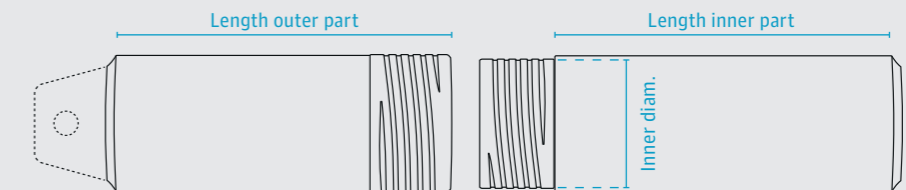

Are you looking for a solid, economical, and versatile packaging tube? You've just found it ...

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.

MicroPackXPress

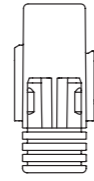

Primary packaging for micro screws

- high quality individual packaging for micro screws offering excellent protection
- linking feature for the flexible composition of sets
- flexible clamping device ensures a firm hold of several diameters of micro screws from 1.5 mm thread diameter
- XPress base design for simple and secure loading/removal of the product
- safe and easy screw removal by direct pick-up
- optimum security and protection for product and end user
- sterile barrier system possible by additional blister



MicroPackXPress — presents your product with the value that a precision medical device deserves. This is ensured by the cap made of clear PETG.

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.



PACKAGING

Base plugs and outer parts are packed separately



MATERIAL

- Outer part: PETG
- Base plugs: Polyethylene (PE)



COLOR

- Outer part in PETG clear
- Base plug in PE natural

USABLE LENGTHS



SHANK DIAM.



PRODUCTION ENVIRONMENT


- Cleanroom production

STANDARD COLOR




Product Range

Article	Suitable for diam.	Usable length fixed
MPXP 8x18-1.5/1	1.5mm	26mm
MPXP 8x18-2/1	1.6 - 2.0mm	26mm
MPXP 8x18-2.35/1	2.1 - 2.35mm	26mm
MPXP 8x18-3/1	2.4 - 3.0mm	26mm

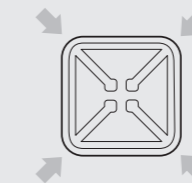

Base plug

PedicleScrewPack

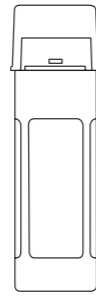

Custom made primary packaging for pedicle screws


- high quality and stable single use package for polyaxial and monoaxial pedicle screws
- developed in cooperation with medical technology companies and clinical specialists
- optimum security and protection for product and end user
- packaging tube for safe removal of the screw using the attaching tool
- fixing of the pedicle screw by product-specific clamping inside the tube
- speedy efficient packaging
- speedy handling
- safe transportation



The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.

Product Range

Article	Description	Usable length	PU
PSP 60-3.2/1 C	short tube with individual clamping	60mm	upon request
PSP 110-3.2/1 C	long tube with individual clamping	110mm	upon request
PSP 27x21/1 B	cap	-	2000

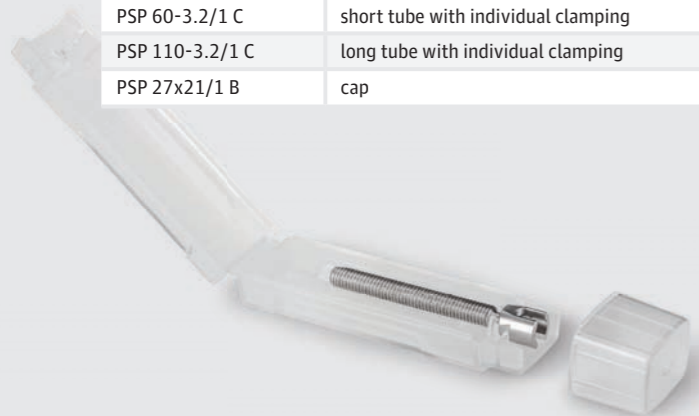

PACKAGING
Tube and cap are packed separately




MATERIAL
— PETG

PRODUCTION ENVIRONMENT
— cleanroom production

STANDARD COLOR
●



Plug Pack

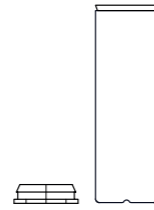

Round packaging tubes with end plug

- stable packaging tube with plug closure
- universal use for products of different dimensions
- square-shaped plug prevents rolling
- security and protection for your product
- speedy efficient packaging
- speedy handling
- safe transportation
- can be used as primary or secondary packaging



PACKAGING


Tubes and plugs are packed separately



MATERIAL

- Tube: Polypropylene (PP)
- Plug: Polyethylene (PE)

USABLE LENGTHS

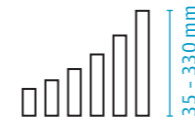

Ø DIAMETER



PRODUCTION ENVIRONMENT

- various sizes available in cleanroom production

STANDARD COLOR




Product Range

Article	Inner diam.	Usable length	PU	Article	Inner diam.	Usable length	PU
TP 7 050 A + P 8	8mm	50mm	3200	TP 20 200 A + P 20	22mm	200mm	900
TP 7 075 A + P 8	8mm	75mm	3200	TP 20 220 A + P 20	22mm	220mm	780
TP 9 040 A + P 9	10mm	35mm	2000	TP 22 050 A + P-P 22	24mm	50mm	1000
TP 9 065 A + P 9	10mm	65mm	2000	TP 22 080 A + P-P 22	24mm	80mm	1300
TP 12 040 A + P 12	14mm	40mm	3000	TP 22 120 A + P-P 22	24mm	120mm	800
TP 12 045 A + P 12	14mm	45mm	3000	TP 22 200 A + P-P 22	24mm	200mm	750
TP 12 050 A + P 12	14mm	50mm	2900	TP 24 035 A + P-P 24	26mm	35mm	1400
TP 12 060 A + P 12	14mm	60mm	2200	TP 24 050 A + P-P 24	26mm	50mm	1900
TP 12 070 A + P 12	14mm	70mm	2200	TP 24 080 A + P-P 24	26mm	80mm	1200
TP 12 080 A + P 12	14mm	80mm	3500	TP 24 120 A + P-P 24	26mm	120mm	800
TP 12 120 A + P 12	14mm	120mm	2000	TP 24 200 A + P-P 24	26mm	200mm	650
TP 12 200 A + P 12	14mm	200mm	2100	TP 26 050 A + P-P 26	28mm	50mm	1700
TP 12 240 A + P 12	14mm	240mm	1600	TP 26 080 A + P-P 26	28mm	80mm	850
TP 14 035 A + P 14	16mm	35mm	2200	TP 26 120 A + P-P 26	28mm	120mm	550
TP 14 050 A + P 14	16mm	50mm	2200	TP 26 200 A + P-P 26	28mm	200mm	500
TP 14 060 A + P 14	16mm	60mm	1800	TP 26 330 A + P-P 26	28mm	330mm	320
TP 14 080 A + P 14	16mm	80mm	2800	TP 28 050 A + P 28	30mm	50mm	1400
TP 14 120 A + P 14	16mm	120mm	1600	TP 28 060 A + P 28	30mm	60mm	1000
TP 14 150 A + P 14	16mm	150mm	1200	TP 28 080 A + P 28	30mm	80mm	800
TP 14 200 A + P 14	16mm	200mm	1500	TP 28 120 A + P 28	30mm	120mm	500
TP 16 035 A + P 16	18mm	35mm	2400	TP 28 200 A + P 28	30mm	200mm	450
TP 16 050 A + P 16	18mm	50mm	1700	TP 30 050 A + P-P 30	33mm	50mm	1200
TP 16 080 A + P 16	18mm	80mm	2000	TP 30 080 A + P-P 30	33mm	80mm	650
TP 16 090 A + P 16	18mm	90mm	1800	TP 30 120 A + P-P 30	33mm	120mm	400
TP 16 120 A + P 16	18mm	120mm	1200	TP 30 150 A + P-P 30	33mm	150mm	350
TP 16 200 A + P 16	18mm	200mm	1200	TP 30 190 A + P-P 30	33mm	190mm	430
TP 16 250 A + P 16	18mm	250mm	900	TP 30 200 A + P-P 30	33mm	200mm	400
TP 18 050 A + P 18	20mm	50mm	1500	TP 30 270 A + P-P 30	33mm	270mm	3000
TP 18 080 A + P 18	20mm	80mm	1700	TP 32 050 A + P-P 32	35mm	50mm	1000
TP 18 120 A + P 18	20mm	120mm	1100	TP 32 080 A + P-P 32	35mm	80mm	650
TP 18 200 A + P 18	20mm	200mm	1000	TP 32 120 A + P-P 32	35mm	120mm	400
TP 18 250 A + P 18	20mm	250mm	800	TP 32 200 A + P-P 32	35mm	200mm	350
TP 20 050 A + P 20	22mm	50mm	1200				
TP 20 060 A + P 20	22mm	60mm	2000				
TP 20 070 A + P 20	22mm	70mm	2000				
TP 20 080 A + P 20	22mm	80mm	1500				
TP 20 120 A + P 20	22mm	120mm	1000				



Oh so simple, oh so practical:
universal packaging tube
with the world's simplest
closure (probably)

QuadroPack

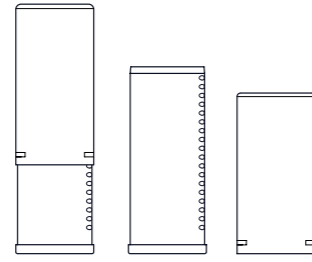

Square telescopic packaging tubes with ratchet mechanism

- infinitely variable telescopic length adjustment
- universal use for products of different dimensions
- security and protection for your product
- speedy efficient packaging
- speedy handling
- safe transportation
- excellent for stacking and storage
- can be used as primary or secondary packaging



PACKAGING
Outer parts and inner parts are packed separately




MATERIAL
— Polypropylene (PP)




COLOR
— PP transparent

USABLE LENGTHS

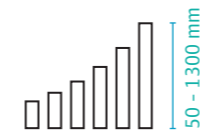

Ø DIAMETER

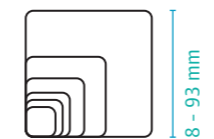

PRODUCTION ENVIRONMENT

— various sizes available in cleanroom production

STANDARD COLOR



Product Range

Article	Inner diam.	Usable length	PU	Article	Inner diam.	Usable length	PU
QP 8 050	8mm	50 - 80mm	2700	QP 26 050	26mm	50 - 60mm	1400
QP 8 080	8mm	80 - 120mm	1700	QP 26 060	26mm	60 - 80mm	1400
QP 8 120	8mm	120 - 200mm	2000	QP 26 080	26mm	80 - 120mm	1000
QP 8 200	8mm	200 - 350mm	2000	QP 26 100	26mm	100 - 160mm	800
QP 10 050	10mm	50 - 80mm	2100	QP 26 120	26mm	120 - 200mm	600
QP 10 080	10mm	80 - 120mm	2600	QP 26 200	26mm	200 - 350mm	600
QP 10 120	10mm	120 - 200mm	1400	QP 26 350	26mm	350 - 620mm	340
QP 10 200	10mm	200 - 350mm	1550	QP 32 050	32mm	50 - 70mm	1400
QP 10 250	10mm	250 - 460mm	1200	QP 32 060	32mm	60 - 80mm	1000
QP 10 300	10mm	300 - 535mm	500	QP 32 080	32mm	80 - 120mm	700
QP 12 050	12mm	50 - 80mm	3000	QP 32 120	32mm	120 - 200mm	450
QP 12 070	12mm	70 - 110mm	2000	QP 32 200	32mm	200 - 350mm	420
QP 12 080	12mm	80 - 120mm	1800	QP 32 350	32mm	350 - 620mm	208
QP 12 120	12mm	120 - 200mm	1000	QP 34 100	34mm	100 - 160mm	450
QP 12 200	12mm	200 - 350mm	1000	QP 35 050	35mm	50 - 70mm	1000
QP 12 250	12mm	250 - 460mm	850	QP 35 060	35mm	60 - 80mm	750
QP 14 050	14mm	50 - 80mm	2300	QP 35 080	35mm	80 - 120mm	600
QP 14 060	14mm	60 - 85mm	2000	QP 35 120	35mm	120 - 200mm	350
QP 14 065	14mm	65 - 100mm	2000	QP 35 200	35mm	200 - 350mm	320
QP 14 080	14mm	80 - 120mm	1300	QP 35 300	35mm	300 - 540mm	210
QP 14 120	14mm	120 - 200mm	800	QP 35 350	35mm	350 - 620mm	168
QP 14 200	14mm	200 - 350mm	800	QP 45 060	45mm	60 - 80mm	450
QP 14 300	14mm	300 - 535mm	500	QP 45 080	45mm	80 - 120mm	300
QP 16 050	16mm	50 - 80mm	2000	QP 45 120	45mm	120 - 200mm	335
QP 16 080	16mm	80 - 120mm	1100	QP 45 200	45mm	200 - 350mm	180
QP 16 120	16mm	120 - 200mm	1400	QP 45 350	45mm	350 - 620mm	99
QP 16 200	16mm	200 - 350mm	1400	QP 45 700	45mm	700 - 1300mm	99
QP 18 050	18mm	50 - 80mm	1500	QP 55 060	55mm	60 - 80mm	300
QP 18 080	18mm	80 - 120mm	1800	QP 55 080	55mm	80 - 120mm	288
QP 18 120	18mm	120 - 200mm	1100	QP 55 120	55mm	120 - 200mm	200
QP 18 200	18mm	200 - 350mm	1100	QP 55 200	55mm	200 - 350mm	120
QP 20 050	20mm	50 - 80mm	1200	QP 55 350	55mm	350 - 620mm	72
QP 20 080	20mm	80 - 120mm	1600	QP 65 080	65mm	80 - 120mm	210
QP 20 120	20mm	120 - 200mm	900	QP 65 120	65mm	120 - 200mm	155
QP 20 200	20mm	200 - 350mm	900	QP 65 200	65mm	200 - 350mm	92
QP 20 350	20mm	350 - 620mm	500	QP 65 300	65mm	300 - 520mm	65
QP 22 050	22mm	50 - 80mm	1200	QP 65 350	65mm	350 - 620mm	56
QP 22 075	22mm	75 - 105mm	1500	QP 78 080	78mm	80 - 120mm	160
QP 22 080	22mm	80 - 120mm	1300	QP 78 120	78mm	120 - 200mm	105
QP 22 120	22mm	120 - 200mm	850	QP-G78150	78mm	150 - 250mm	80
QP 22 200	22mm	200 - 350mm	840	QP 78 220	78mm	220 - 350mm	56
QP 24 050	24mm	50 - 80mm	850	QP 78 320	78mm	320 - 450mm	42
QP 24 080	24mm	80 - 120mm	1000	QP 93 080	93mm	80 - 120mm	90
QP 24 120	24mm	120 - 200mm	650	QP 93 120	93mm	120 - 200mm	75
QP 24 200	24mm	200 - 350mm	600	QP 93 200	93mm	200 - 350mm	48
				QP 93 280	93mm	280 - 450mm	30

QP-G 78 150 with recessed bottom



Square, practical, easily closable – thanks to ratchet length adjustment. It's no accident that our cuboid multi-talent is among the world's best-selling packaging tubes

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.

QuadroPack TipLine



Protective packaging tubes for surgical instruments and applications

- square protective packaging tube with tapered ends
- optimum protection for long surgical instruments and various other applications
- tapered ends reduce the risk of the content breaking through
- universal use for products of different dimensions
- security and protection for your product
- speedy and efficient packing of products
- safe transportation
- excellent for stacking and storage
- can be used as primary or secondary packaging



PACKAGING
Outer parts and inner parts are packed separately




MATERIAL
- PETG




COLOR
- PETG

USABLE LENGTHS

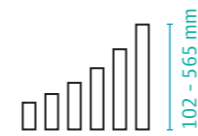

Ø DIAMETER

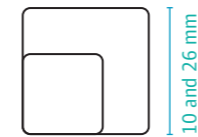

PRODUCTION ENVIRONMENT

- cleanroom production


STANDARD COLOR



Product Range


Article	Ø	Usable length	PU
QPTL 10 050 A + 10 080 I	10mm	102-132mm	2000
QPTL 10 050 A + 10 120 I	10mm	142-172mm	2000
QPTL 10 050 A + 10 200 I	10mm	222-252mm	2000
QPTL 10 050 A + 10 250 I	10mm	272-302mm	2000
QPTL 10 050 A + 10 300 I	10mm	322-352mm	2000
QPTL 10 200 A + 10 080 I	10mm	302mm	1000
QPTL 10 200 A + 10 120 I	10mm	342mm	1000
QPTL 10 200 A + 10 200 I	10mm	422mm	1000
QPTL 10 200 A + 10 250 I	10mm	472mm	1000
QPTL 10 200 A + 10 300 I	10mm	522mm	1000
QPTL 26 050 A + 26 120 I	26mm	161-185mm	1000
QPTL 26 050 A + 26 200 I	26mm	241-265mm	1000
QPTL 26 050 A + 26 350 I	26mm	391-415mm	1000
QPTL 26 050 A + 26 500 I	26mm	541-565mm	1000

The usable length describes the maximum internal length of the packaging.

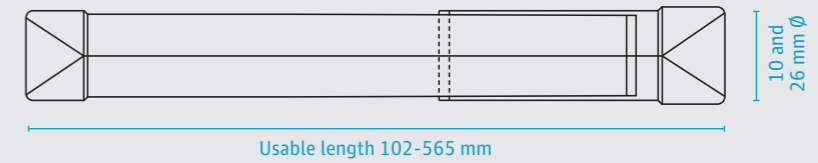

All the advantages of our QuadroPack – plus optimized protection for long surgical instruments and other applications

Rectangular Pack

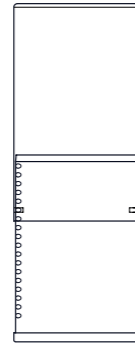

Rectangular telescopic packaging tubes with ratchet mechanism

- infinitely variable telescopic length adjustment
- universal use for products of different dimensions
- security and protection for your product
- speedy efficient packaging
- speedy handling
- safe transportation
- excellent for stacking and storage
- can be used as primary or secondary packaging



PACKAGING
Outer parts and inner parts are packed separately




MATERIAL
- Polyethylene (PE)
- Polypropylene (PP)




COLOR
- PE natural
- PP transparent

USABLE LENGTHS

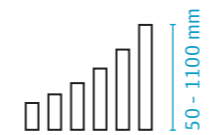

DIMENSION

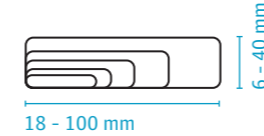

PRODUCTION ENVIRONMENT

- various sizes available in cleanroom production

STANDARD COLOR



Product Range

Article	Inner depth	Inner width	Usable length adjustable	PU	Article	Inner depth	Inner width	Usable length adjustable	PU
RP 6x20/050	6mm	20mm	50-80mm	2000	RP 19x26/245	19mm	26mm	245-430mm	550
RP 6x20/080	6mm	20mm	80-120mm	1000	RP 19x31/270	19mm	31mm	270-470mm	400
RP 6x20/120	6mm	20mm	120-200mm	1000	RP 19x31/295	19mm	31mm	295-520mm	360
RP 6x20/200	6mm	20mm	200-360mm	1000	RP 25x33/080	25mm	33mm	80-120mm	700
RP 8x36/050	8mm	36mm	50-80mm	1600	RP 26x31/270	26mm	31mm	270-470mm	330
RP 8x36/080	8mm	36mm	80-120mm	1000	RP 26x42/060	26mm	42mm	60-80mm	700
RP 8x36/138	8mm	36mm	138-240mm	450	RP 26x42/080	26mm	42mm	80-120mm	-
RP 8x36/240	8mm	36mm	240-420mm	480	RP 26x42/120	26mm	42mm	120-200mm	400
RP 8x36/420	8mm	36mm	420-780mm	190	RP 26x42/220	26mm	42mm	220-380mm	320
RP 8x36/820	8mm	36mm	820-1100mm	190	RP 26x42/335	26mm	42mm	335-580mm	200
RP 9x45/080	9mm	45mm	80-120mm	550	RP 29x36/080	29mm	36mm	80-120mm	700
RP 9x45/120 I	9mm	45mm	120-160mm	750	RP 30x65/055	30mm	65mm	55-70mm	650
RP 9x45/160 I	9mm	45mm	160-200mm	270	RP 30x65/070	30mm	65mm	70-100mm	520
RP 10x18/080	10mm	18mm	80-120mm	1200	RP 30x65/100	30mm	65mm	100-150mm	350
RP 10x18/120	10mm	18mm	120-200mm	1500	RP 30x65/130	30mm	65mm	130-200mm	385
RP 10x27/070 I	10mm	27mm	70-110mm	850	RP 30x65/160	30mm	65mm	160-260mm	300
RP 10x27/080	10mm	27mm	80-120mm	850	RP 30x65/250	30mm	65mm	250-450mm	180
RP 10x33/120	10mm	33mm	120-200mm	500	RP 30x65/300	30mm	65mm	300-550mm	160
RP 10x33/320	10mm	33mm	320-550mm	240	RP 31x56/445	31mm	56mm	445-780mm	100
RP 12x26/200	12mm	26mm	200-360mm	480	RP 33x42/050	33mm	42mm	50-80mm	800
RP 12x36/360	12mm	36mm	360-620mm	330	RP 33x42/335	33mm	42mm	335-580mm	165
RP 14x36/080	14mm	36mm	80-120mm	650	RP 36x62/445	36mm	62mm	445-780mm	80
RP 16x36/050	16mm	36mm	50-80mm	700	RP 40x65/310	40mm	65mm	310-310mm	105
RP 16x36/080	16mm	36mm	80-120mm	600	RP 40x65/540	40mm	65mm	540-540mm	110
RP 16x36/120	16mm	36mm	120-200mm	350	RP 40x75/055	40mm	75mm	55-70mm	450
RP 16x36/200	16mm	36mm	200-360mm	350	RP 40x75/070	40mm	75mm	70-100mm	350
RP 16x36/360 I	16mm	36mm	360-360mm	200	RP 40x75/100	40mm	75mm	100-150mm	320
RP 16x55/120	16mm	55mm	120-200mm	450	RP 40x75/150	40mm	75mm	150-250mm	200
RP 16x55/200	16mm	55mm	200-360mm	370	RP 40x75/210	40mm	75mm	210-360mm	150
RP 18x42/050	18mm	42mm	50-80mm	1000	RP 40x75/250	40mm	75mm	250-450mm	120
RP 18x42/080	18mm	42mm	80-120mm	850	RP 40x75/400	40mm	75mm	400-700mm	70
RP 18x42/120	18mm	42mm	120-200mm	450	RP 40x100/150	40mm	100mm	150-250mm	132
RP 18x42/200	18mm	42mm	200-360mm	408					
RP 18x42/240	18mm	42mm	240-420mm	180					
RP 18x42/335	18mm	42mm	335-580mm	280					
RP 18x42/620	18mm	42mm	620-1100mm	180					




Not everything is round or square – which is why we created this practical packaging tube for everything that's on the flatter side

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.


Perfect protection for long products

RectangularPack TipLine

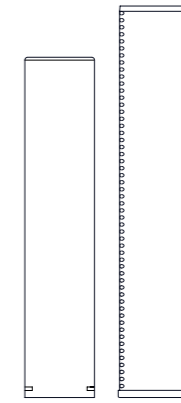

Rectangular telescopic packaging tubes with tapered ends

- universal use for long products
- fixed usable length due to length stop
- security and protection for your product
- speedy and efficient packing of products
- safe transportation
- excellent for stacking and storage
- can be used as primary or secondary packaging

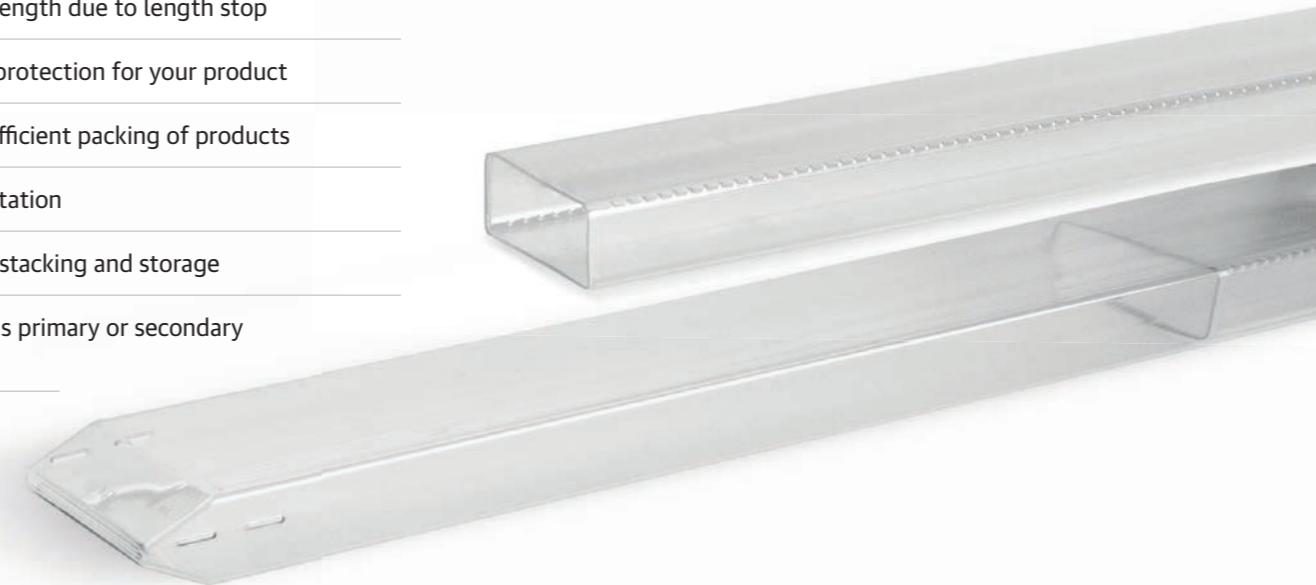

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.

Product Range

Inner part (I)	Outer part (A)	Ø	Usable length	PU (I)	PU (A)
RPTL 18x46/450 I	RPTL 18x46/382 A	18x46mm	869mm	upon request	100

The usable length describes the maximum internal length of the packaging.



PACKAGING
Outer parts and inner parts are packed separately




MATERIAL
- PETG

PRODUCTION ENVIRONMENT
- cleanroom production

STANDARD COLOR


TelePack

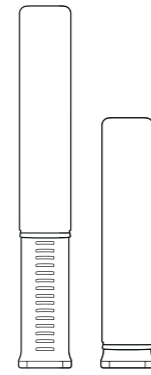

Round telescopic packaging tubes with ratchet length adjustment

- telescopic length extension, ratchet mechanism locks package into extended position
- square end prevents rolling
- universal use for products of different dimensions
- security and protection for your product
- speedy efficient packaging
- speedy handling
- safe transportation
- can be used as primary or secondary packaging



PACKAGING
Outer parts and inner parts are packed separately




MATERIAL
— Polypropylene (PP)



COLOR
— PP transparent

USABLE LENGTH

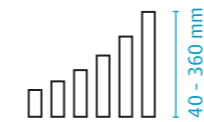

Ø DIAMETER



PRODUCTION ENVIRONMENT

— various sizes available in cleanroom production

STANDARD COLOR



Product Range

Article	Inner diam.	Usable length (adjustable)	PU
TP 7 050	7mm	50 - 90mm	3200
TP 7 075	7mm	75 - 115mm	3200
TP 9 040	9mm	40 - 65mm	2000
TP 9 065	9mm	65 - 110mm	2000
TP 12 050	12mm	50 - 80mm	2900
TP 12 080	12mm	80 - 120mm	3500
TP 12 120	12mm	120 - 200mm	2000
TP 12 200	12mm	200 - 360mm	2100
TP 14 050	14mm	50 - 80mm	2200
TP 14 080	14mm	80 - 120mm	2800
TP 14 120	14mm	120 - 200mm	1600
TP 14 200	14mm	200 - 360mm	1500
TP 16 050	16mm	50 - 80mm	1700
TP 16 080	16mm	80 - 120mm	2000
TP 16 120	16mm	120 - 200mm	1200
TP 16 200	16mm	200 - 360mm	1200
TP 18 050	18mm	50 - 80mm	1500
TP 18 080	18mm	80 - 120mm	1700
TP 18 120	18mm	120 - 200mm	1100
TP 18 200	18mm	200 - 360mm	1000
TP 20 050	20mm	50 - 80mm	1200
TP 20 080	20mm	80 - 120mm	1500
TP 20 120	20mm	120 - 200mm	1000
TP 20 200	20mm	200 - 360mm	900
TP 22 050	22mm	50 - 80mm	1000
TP 22 080	22mm	80 - 120mm	1300
TP 22 120	22mm	120 - 200mm	800
TP 22 200	22mm	200 - 360mm	750
TP 24 050	24mm	50 - 80mm	1900
TP 24 080	24mm	80 - 120mm	1200
TP 24 120	24mm	120 - 200mm	800
TP 24 200	24mm	200 - 360mm	650
TP 26 050	26mm	50 - 80mm	1700
TP 26 080	26mm	80 - 120mm	850
TP 26 120	26mm	120 - 200mm	550
TP 26 200	26mm	200 - 360mm	500
TP 28 050	28mm	50 - 80mm	1400
TP 28 080	28mm	80 - 120mm	800
TP 28 120	28mm	120 - 200mm	500
TP 28 200	28mm	200 - 360mm	450
TP 30 050	30mm	50 - 80mm	1200
TP 30 080	30mm	80 - 120mm	650
TP 30 120	30mm	120 - 200mm	400
TP 30 200	30mm	200 - 360mm	400
TP 32 050	32mm	50 - 80mm	1000
TP 32 080	32mm	80 - 120mm	650
TP 32 120	32mm	120 - 200mm	400
TP 32 200	32mm	200 - 360mm	350



The telescopic length adjustment was the inspiration for the name of this practical and versatile packaging tube – a sought-after classic in our product range

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.

TipTube

Adjustable length package for bone nails with a bayonet locking feature

- high-quality, solid individual packaging for bone nails of various geometries (e.g. femoral, humeral, tibial)
- optimum security and protection for product and user
- secure positioning of bone nail
- minimal contact of packaging and bone nail
- easy handling due to bayonet lock
- secure handover of the packaged product
- can be used as primary or secondary packaging




PACKAGING
Outer parts and inner parts are packed separately




MATERIAL
— PETG

USABLE LENGTHS

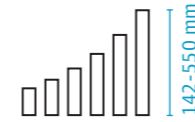

Ø DIAMETER

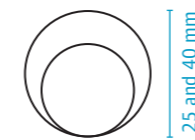

PRODUCTION ENVIRONMENT

— cleanroom production


STANDARD COLOR



Product Range


Article	Inner diam.	Usable length	PU
TI 25 110	25mm	142-202mm	upon request
TI 40 280	40mm	330-550mm	upon request

The usable length describes the maximum internal length of the packaging.

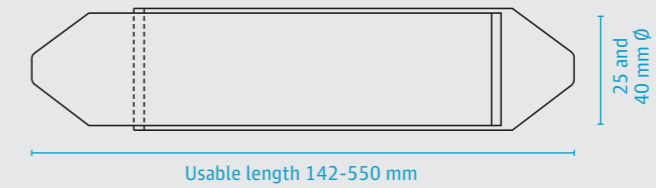

This length-adjustable bone nail packaging provides optimum security and protection for product and user

TopPack XPress

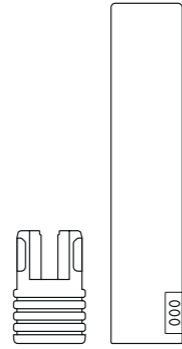

Single packaging with short shankholding base

- XPress base design for simple and secure loading/removal of the product
- internal ribbing molded in plug keeps the tool perfectly positioned
- user-friendly function ensures easier and safer handling of sensitive and sharp-edged products
- security and protection for your product
- speedy efficient packaging
- speedy handling
- safe transportation
- can be used as primary packaging



 A secure grip with an attractive design

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.



PACKAGING
Base plugs and outer parts are packed separately




MATERIAL
- Outer part: Polypropylene (PP)
- Base plugs: Polyethylene (PE)



COLOR
- Outer part in PP transparent
- Base plug in PE natural

USABLE LENGTHS

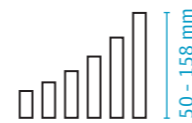

SHANK DIAM.



PRODUCTION ENVIRONMENT

- various sizes available in cleanroom production

STANDARD COLOR




Product Range

Base plug

Article	Part type	Diam.	Shank diam.	Shank diam. (in.)	Inner length	PU
TPKXP 8x18-2/1 I	base plug	8mm	2mm		18mm	1000
TPKXP 8x18-2.35/1 I	base plug	8mm	2.35mm		18mm	1000
TPKXP 8x18-3/1 I	base plug	8mm	3mm		18mm	1000
TPKXP 8x18-1/8"/1 I	base plug	8mm	3.175mm	1/8"	18mm	1000
TPKXP 8x18-4/1 I	base plug	8mm	4mm		18mm	1000
TPKXP 8x18-3/16"/1 I	base plug	8mm	4.762mm	3/16"	18mm	1000
TPKXP 8x18-5/1 I	base plug	8mm	5mm		18mm	1000
TPKXP 8x18-6/1 I	base plug	8mm	6mm		18mm	1000
TPKXP 8x18-1/4"/1 I	base plug	8mm	6.35mm	1/4"	18mm	1000
TPKXP 8x18-7/1 I	base plug	8mm	7mm		18mm	1000
TPKXP 8x18-8/1 I	base plug	8mm	8mm		18mm	1000
TPKXP 13x25-2.35/1 I	base plug	13mm	2.35mm		25mm	1000
TPKXP 13x25-3/1 I	base plug	13mm	3mm		25mm	1000
TPKXP 13x25-4/1 I	base plug	13mm	4mm		25mm	1000
TPKXP 13x25-5/1 I	base plug	13mm	5mm		25mm	1000
TPKXP 13x25-6/1 I	base plug	13mm	6mm		25mm	1000
TPKXP 13x25-1/4"/1 I	base plug	13mm	6.35mm	1/4"	25mm	1000
TPKXP 13x25-7/1 I	base plug	13mm	7mm		25mm	1000
TPKXP 13x25-8/1 I	base plug	13mm	8mm		25mm	1000
TPKXP 13x25-9/1 I	base plug	13mm	9mm		25mm	1000
TPKXP 13x25-3/8"/1 I	base plug	13mm	9.525mm	3/8"	25mm	1000
TPKXP 13x25-10/1 I	base plug	13mm	10mm		25mm	1000
TPKXP 13x25-7/16"/1 I	base plug	13mm	11.11mm	7/16"	25mm	1000
TPKXP 13x25-12/1 I	base plug	13mm	12mm		25mm	1000
TPKXP 13x25-1/2"/1 I	base plug	13mm	12.7mm	1/2"	25mm	1000
TPKXP 18x34-9/1 I	base plug	18mm	9mm		34mm	250
TPKXP 18x34-12/1 I	base plug	18mm	12mm		34mm	250
TPKXP 18x34-13/1 I	base plug	18mm	13mm		34mm	250
TPKXP 18x34-14/1 I	base plug	18mm	14mm		34mm	250
TPKXP 18x34-15/1 I	base plug	18mm	15mm		34mm	250
TPKXP 18x34-9/16"/1 I	base plug	18mm	14.287mm	9/16"	34mm	250
TPKXP 18x34-5/8"/1 I	base plug	18mm	15.875mm	5/8"	34mm	250
TPKXP 18x34-16/1 I	base plug	18mm	16mm		34mm	250
TPKXP 20x40-18/1 I	base plug	20mm	18mm		40mm	250
TPKXP 20x40-3/4"/1 I	base plug	20mm	19.05mm	3/4"	40mm	250
TPKXP 20x40-20/1 I	base plug	20mm	20mm		40mm	250

Outer part see next page ↗

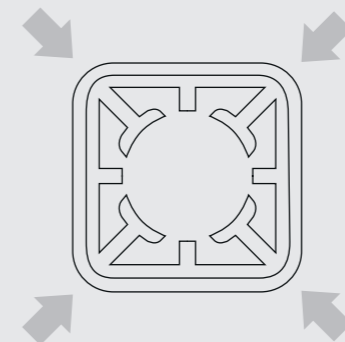

Base plug



PACKAGING

Base plugs and outer parts are packed separately



MATERIAL


- Outer part: Polypropylene (PP)
- Base plugs: Polyethylene (PE)




COLOR

- Outer part in PP transparent
- Base plug in PE natural

USABLE LENGTHS

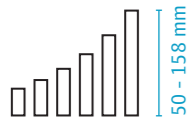

SHANK DIAM.

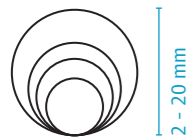

PRODUCTION ENVIRONMENT

- various sizes available in cleanroom production

STANDARD COLOR




Product Range

Outer part

Article	Part type	Diameter	Inner length	Suitable base plug	Usable length	PU
TPK 8x42 A	outer part	8mm	42mm	TPKXP 8	50	1000
TPK 8x52 A	outer part	8mm	52mm	TPKXP 8	60	1000
TPK 8x70 A	outer part	8mm	70mm	TPKXP 8	77	1000
TPK 8x85/2 A	outer part	8mm	85mm	TPKXP 8	92	2300
TPK 8x97 A	outer part	8mm	97mm	TPKXP 8	105	1000
TPK 8x97/2 A	outer part	8mm	97mm	TPKXP 8	105	1000
TPK 13x56 A	outer part	13mm	56mm	TPKXP 13	66	500
TPK 13x74 A	outer part	13mm	74mm	TPKXP 13	84	500
TPK 13x88 A	outer part	13mm	88mm	TPKXP 13	98	500
TPK 13x97 A	outer part	13mm	97mm	TPKXP 13	108	500
TPK 13x110 A	outer part	13mm	110mm	TPKXP 13	120	500
TPK 18x83 A	outer part	18mm	83mm	TPKXP 18	98	500
TPK 18x113 A	outer part	18mm	113mm	TPKXP 18	129	500
TPK 18x130 A	outer part	18mm	130mm	TPKXP 18	144	250
TPK 18x145/2 A	outer part	18mm	145mm	TPKXP 18	158	500
TPK 20x89 A	outer part	20mm	89mm	TPKXP 20	108	250
TPK 20x130 A	outer part	20mm	130mm	TPKXP 20	149	250



TubeInTube®

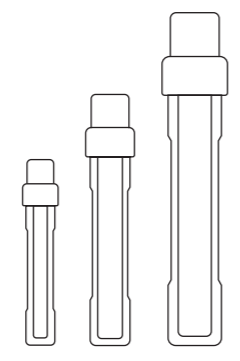

PRODUCT UNDER DEVELOPMENT

Double sterile barrier packaging

- double sterile barrier packaging system
- outer and inner tube with sterile barrier and tamper proof closure
- creates a secure, effective sterile barrier
- protects and safeguards your implant's functionality
- makes sterilization in the hospital obsolete
- optimizes the process and simplifies handling in operating rooms
- saves storage space



The new TubeInTube® double sterile barrier packaging system makes operating room sterilization obsolete

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.



MATERIAL


- Inner and outer tube: PETG
- Inner and outer plug: combination of PE and TPE
- Tamper proof closure: PE



COLOR

- PETG clear

USABLE LENGTHS

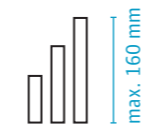



PRODUCTION ENVIRONMENT

- cleanroom production

Product Range


Article	Inner tube Ø	Outer tube Ø	Overall length	Screw thread Ø min.-max.	Screw length max.
TT 10 040 K*	10	16	97	2.0, 2.5, 3.0, 3.5	40
TT 10 080 K*	10	16	137	4.0, 5.0, 6.0	80
TT 10 160 K*	10	16	217	7.0, 8.0	160

*planned dimensions; product in development stage

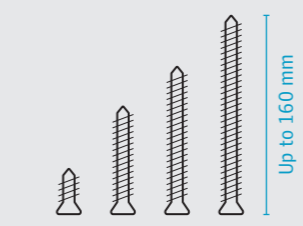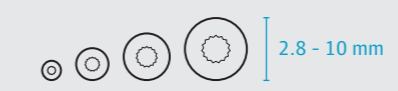


BONE SCREW SIZES

TubeInTube® is suitable for bone screws with thread diameters from 2 to 8 mm (screw head diameters from 2.8 to 10 mm) and in lengths up to 160 mm.

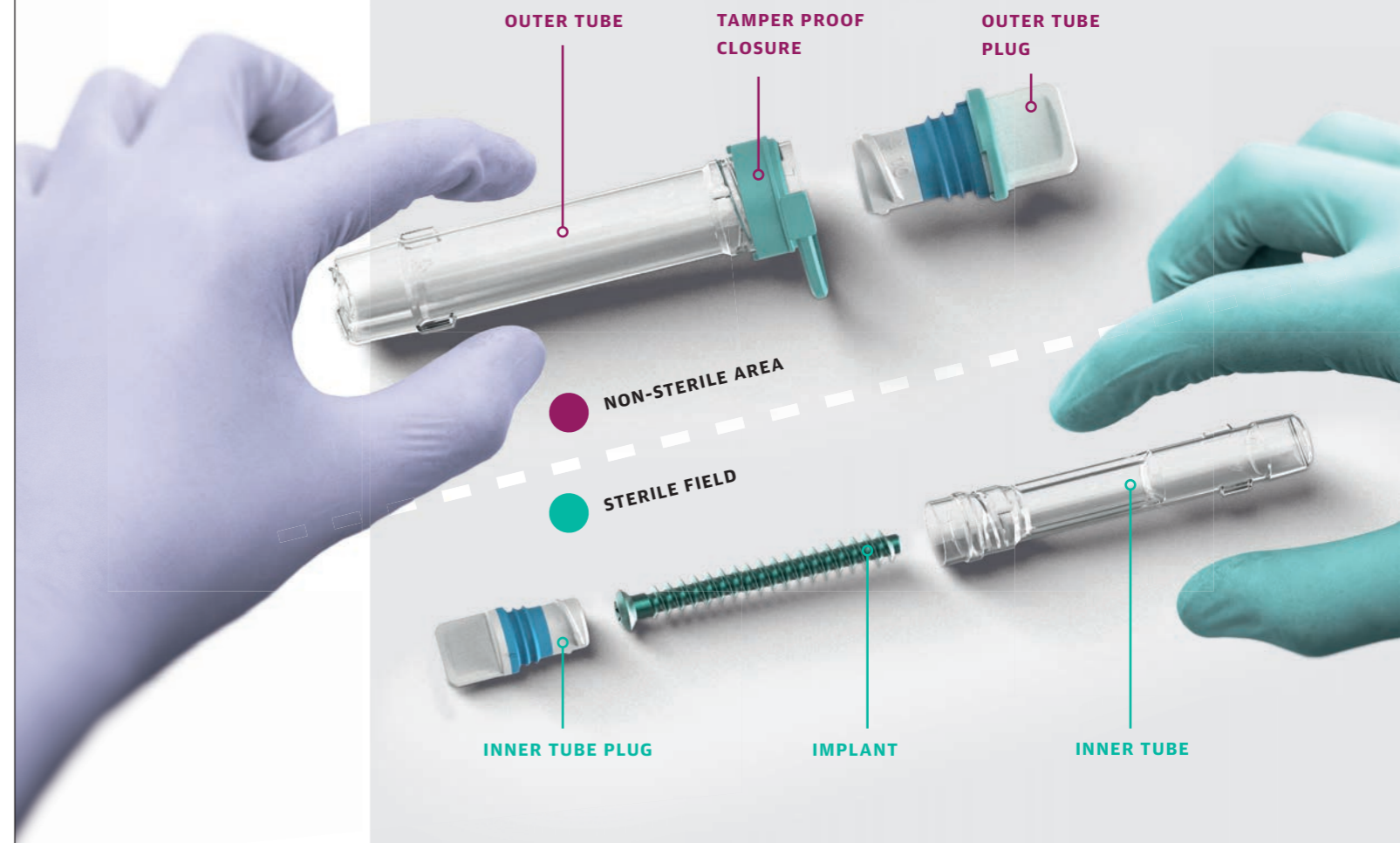

TwistPackPlus

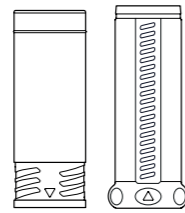

More safety. More efficiency. The new generation of TwistPack packaging tubes convince with user-friendly twist lock closure

- maximum flexibility through universal length adjustment
- optimal product protection
- easy opening and closing through telescopic adjustment
- firm hold force through our click groove system
- hexagon tube bottom prevents rolling
- speedy handling
- safe transportation
- durable and reusable
- can be used as primary or secondary packaging



Need packaging for elongated or cylindrical products? Meet the solution. Simply adjust the length, turn – et voila!

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.



PACKAGING

Outer parts and inner parts are packed separately:
TX 85 and TX 105 outer parts and inner parts are packed assembled



MATERIAL


- Polypropylene (PP)
- Polyethylene (PE)



COLOR

- PP transparent
- PE natural

USABLE LENGTHS

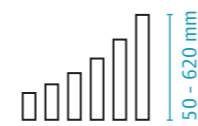

Ø DIAMETER



PRODUCTION ENVIRONMENT

- various sizes available in cleanroom production

STANDARD COLOR




Product Range

Article	Inner diam.	Usable length	PU	Article	Inner diam.	Usable length	PU
TX 10 050	10mm	50 - 70mm	2000	TX 35 050	35mm	50 - 60mm	450
TX 10 080	10mm	80 - 120mm	1700	TX 35 060	35mm	60 - 80mm	600
TX 10 120	10mm	120 - 200mm	1000	TX 35 080	35mm	80 - 120mm	550
TX 13 050	13mm	50 - 70mm	2100	TX 35 120	35mm	120 - 200mm	350
TX 13 080	13mm	80 - 120mm	1300	TX 35 200	35mm	200 - 350mm	250
TX 13 120	13mm	120 - 200mm	1700	TX 35 350	35mm	350 - 620mm	150
TX 13 200	13mm	200 - 350mm	1350	TX 35 500*	35mm	500 - 750mm	250
TX 16 050	16mm	50 - 70mm	1500	TX 35 700*	35mm	700 - 1050mm	150
TX 16 080	16mm	80 - 120mm	2000	TX 40 050	40mm	50 - 60mm	700
TX 16 120	16mm	120 - 200mm	1200	TX 40 060	40mm	60 - 80mm	600
TX 16 200	16mm	200 - 350mm	1000	TX 40 080	40mm	80 - 120mm	450
TX 16 350	16mm	350 - 620mm	500	TX 40 120	40mm	120 - 200mm	270
TX 16 410	16mm	410 - 770mm	500	TX 40 200	40mm	200 - 350mm	208
TX 20 050	20mm	50 - 70mm	950	TX 40 350	40mm	350 - 620mm	120
TX 20 080	20mm	80 - 120mm	1500	TX 45 050	45mm	50 - 60mm	280
TX 20 120	20mm	120 - 200mm	900	TX 45 060	45mm	60 - 80mm	500
TX 20 200	20mm	200 - 350mm	700	TX 45 080	45mm	80 - 120mm	300
TX 22 050	22mm	50 - 60mm	1000	TX 45 120	45mm	120 - 200mm	250
TX 22 080	22mm	80 - 120mm	1300	TX 45 200	45mm	200 - 350mm	180
TX 22 120	22mm	120 - 200mm	650	TX 45 350	45mm	350 - 620mm	100
TX 22 200	22mm	200 - 350mm	600	TX 45 500*	45mm	500 - 750mm	180
TX 22 350	22mm	350 - 620mm	390	TX 45 700*	45mm	700 - 1050mm	100
TX 22 410	22mm	410 - 770mm	300	TX 50 050	50mm	50 - 60mm	400
TX 22 500*	22mm	500 - 750mm	600	TX 50 060	50mm	60 - 80mm	350
TX 22 700*	22mm	700 - 1050mm	390	TX 50 080	50mm	80 - 120mm	300
TX 24 050	24mm	50 - 60mm	810	TX 50 120	50mm	120 - 200mm	190
TX 24 080	24mm	80 - 120mm	1000	TX 50 200	50mm	200 - 350mm	140
TX 24 120	24mm	120 - 200mm	600	TX 50 350	50mm	350 - 620mm	80
TX 24 200	24mm	200 - 350mm	500	TX 55 050	55mm	50mm	400
TX 26 050	26mm	50 - 60mm	800	TX 55 060	55mm	60 - 80mm	200
TX 26 080	26mm	80 - 120mm	900	TX 55 080	55mm	80 - 120mm	200
TX 26 120	26mm	120 - 200mm	550	TX 55 120	55mm	120 - 200mm	160
TX 26 200	26mm	200 - 350mm	450	TX 55 200	55mm	200 - 350mm	110
TX 26 350	26mm	350 - 620mm	260	TX 55 350	55mm	350 - 620mm	68
TX 26 410	26mm	410 - 770mm	200	TX 65 050	65mm	50mm	220
TX 26 500*	26mm	500 - 750mm	450	TX 65 060	65mm	60 - 80mm	220
TX 26 700*	26mm	700 - 1050mm	260	TX 65 080	65mm	80 - 120mm	170
TX 28 050	28mm	50 - 60mm	600	TX 65 120	65mm	120 - 200mm	128
TX 28 080	28mm	80 - 120mm	900	TX 65 200	65mm	200 - 350mm	86
TX 28 120	28mm	120 - 200mm	550	TX 65 350	65mm	350 - 620mm	46
TX 28 200	28mm	200 - 350mm	405	TX 85 080	85mm	80 - 100mm	100
TX 30 050	30mm	50 - 60mm	1000	TX 85 120	85mm	120 - 180mm	65
TX 30 080	30mm	80 - 120mm	750	TX 85 200	85mm	200 - 350mm	36
TX 30 120	30mm	120 - 200mm	450	TX 85 350	85mm	350 - 550mm	20
TX 30 200	30mm	200 - 350mm	320	TX 85 550	85mm	550 - 1000mm	18
TX 30 350	30mm	350 - 620mm	200	TX 105 080	105mm	80 - 100mm	72
TX 30 410	30mm	410 - 770mm	160	TX 105 120	105mm	120 - 180mm	48
				TX 105 200	105mm	200 - 350mm	30

*assembled from three parts (two outer parts and one middle part)

Consumer Box

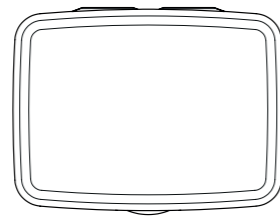

Plastic boxes with hinged lid and integrated latch

- high quality image
- universal application
- polished lid surface provides high transparency for product visibility
- latch integrated into box contours
- durable hinges
- excellent surface for printing and labeling



MATERIAL
- Polypropylene (PP)




COLOR
- PP transparent

HEIGHTS



DIMENSION

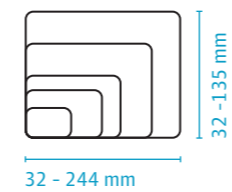

PRODUCTION ENVIRONMENT

- cleanroom production available upon request

STANDARD COLOR



Product Range


Article	Version	Description	Inner length	Inner width	Inner height	PU
CB 34x34x8/1	1		32mm	32mm	8mm	1000
CB 56x42x10/1	1		55mm	39mm	10mm	400
CB 68x50x13/1	1		66mm	47mm	13mm	470
CB 95x70x17/1	1		92mm	65mm	17mm	340
CB 95x70x17/3	3	w/ flat lid	92mm	65mm	17mm	340
CB 115x85x26/1	1		111mm	78mm	26mm	172
CB 138x102x30/1	1		134mm	95mm	30mm	88
CB 150x50x30/1	1		147mm	44mm	29mm	180
CB 180x145x45/1	1		173mm	135mm	45mm	80
CB 246x88x20/1	1		244mm	83mm	16mm	100



Simply brilliant, brilliantly simple: This high-quality plastic box with integrated latch and polished lid surface is the perfect sales packaging

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.

SplitPack

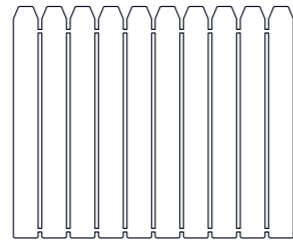

Multi-compartment packaging system with individual detachable units

- versatile packaging combining single unit with multiple packaging
- SplitPack: 2-part version with separate base and lid (SP 4x40 and SP 4x60 are one-piece version)
- SplitPack with shankholding: with integrated shankholding feature in lid
- detachable units make this product very efficient and user-friendly
- no packing and repacking of single items necessary
- each unit is reusable with hinged lid and snap-closure
- speedy handling
- excellent for stacking and storage



PACKAGING
Outer parts and inner parts are packed separately




MATERIAL
— Polypropylene (PP)



COLOR
— PP transparent

USABLE LENGTHS

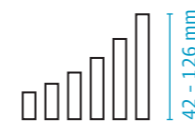

Ø DIAMETER



PRODUCTION ENVIRONMENT

— cleanroom production available upon request

STANDARD COLOR



Product Range

SplitPack

Article	Number of units	Inner diam. unit	Inner length	PU
SP 4x40/2	10	4mm	42mm	700
SP 6x40/2	10	6mm	42mm	520
SP 6x50/2	10	6mm	50mm	780
SP 6x66/2	10	6mm	66mm	625
SP 6x80/2	10	6mm	80mm	500
SP 6x100/2	10	6mm	103mm	405
SP 8x66/2	10	8mm	67mm	600
SP 8x80/2	10	8mm	82mm	500
SP 8x100/2	10	8mm	104mm	420
SP 8x125/2	10	8mm	126mm	340
SP 10x66/2	10	10mm	67mm	420
SP 10x80/2	10	10mm	82mm	360
SP 10x100/2	10	10mm	104mm	300
SP 10x125/2	10	10mm	126mm	250
SP 12x80/2	5	12mm	85mm	525
SP 12x110/2	5	12mm	114mm	400
SP 16x80/2	5	16mm	84mm	330
SP 16x110/2	5	16mm	114mm	250


SplitPack with shankholding

Article	Number of units	Inner diam. unit	Shankholding diam	Inner length	PU
SP 6x50-2/2	10	6mm	2mm	50mm	780
SP 6x50-2.35/2	10	6mm	2.35mm	50mm	780
SP 6x50-2.8/2	10	6mm	2.35mm	50mm	780
SP 6x50-3/2	10	6mm	3mm	50mm	780
SP 6x66-2.35/2	10	6mm	2.35mm	66mm	625
SP 6x66-3/2	10	6mm	3mm	66mm	625
SP 8x66-3/2	10	8mm	3mm	67mm	600
SP 10x80-4/2	10	10mm	4mm	82mm	360
SP 10x80-6/2	10	10mm	6mm	82mm	360
SP 16x80-6/2	5	16mm	6mm	84mm	330



The only packaging we sell that's allowed to be broken. For safe transport and protection of multiple shank tools

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.

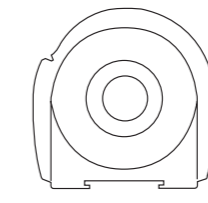

TapePack

Functional single packaging for kinesiology tapes with roll-off function

- high-quality, stable packaging especially for kinesiology tapes
- suitable for space-saving wall assembly by simply hanging onto a rail
- with roll-off function
- recessed grip on the one side, round label surface on the other side
- closure of the two half-shells by snap mechanism
- for tape rolls up to max. 70 mm diameter and 50 mm width


ROLL IT.
TAKE IT.
TAPE IT.



The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.

Lieferprogramm

Article	Dimension	Description
TapePack	86x83x63mm	for tape rolls up to max. 70 mm diameter, 50 mm width

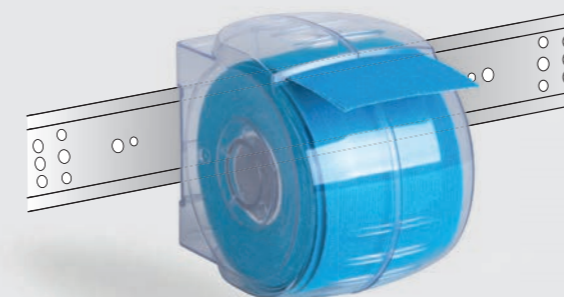

PACKAGING
half-shells are packed separately




MATERIAL
- ABS

PRODUCTION ENVIRONMENT
- Industrial standard

STANDARD COLOR



UniBox

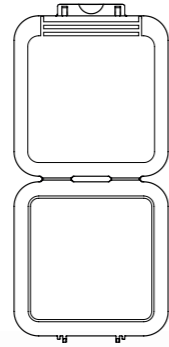

Plastic storage boxes with hinged lid and secure closure

- with secure press closure and hinged lid
- speedy efficient packaging
- speedy handling
- safe transportation
- excellent for stacking and storage
- durable and reusable
- can be used as primary or secondary packaging



MATERIAL
- Polypropylene (PP)




COLOR
- PP transparent

HEIGHTS



DIMENSIONS

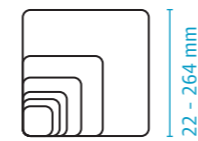

PRODUCTION ENVIRONMENT

- cleanroom production available upon request

STANDARD COLOR



Product Range

UniBox square

Article	Length	Width	Height	PU
UB 22x8/1	22mm	22mm	8mm	1500
UB 35x4/1	37mm	37mm	4mm	1000
UB 35x8/1	36mm	36mm	8mm	1000
UB 35x12/1	35mm	35mm	12mm	1000
UB 35x18/1	35mm	35mm	18mm	750
UB 35x25/1	35mm	35mm	25mm	600
UB 35x40/1	35mm	35mm	40mm	400
UB 47x6/1	49mm	49mm	6mm	750
UB 47x19/1	47mm	47mm	19mm	730
UB 47x40/1	43mm	43mm	40mm	320
UB 58x6/1	63mm	63mm	6mm	450
UB 58x12/1	62mm	62mm	12mm	300
UB 58x25/1	60mm	60mm	25mm	400
UB 58x40/1	58mm	58mm	40mm	330
UB 67x6/1	70mm	70mm	6mm	350
UB 67x12/1	69mm	69mm	12mm	250
UB 67x26/1	67mm	67mm	26mm	300
UB 67x40/1	67mm	67mm	40mm	225
UB 85x6/1	92mm	92mm	6mm	400
UB 85x12/1	91mm	91mm	12mm	250
UB 85x25/1	88mm	88mm	25mm	150
UB 85x40/1	86mm	86mm	40mm	150
UB 100x6/1	105mm	105mm	6mm	300
UB 100x12/1	103mm	103mm	12mm	200
UB 100x25/1	101mm	101mm	25mm	250
UB 100x40/1	100mm	100mm	40mm	110
UB 125x6/1	133mm	133mm	6mm	180
UB 125x12/1	132mm	132mm	12mm	120
UB 125x25/1	130mm	130mm	25mm	150
UB 125x40/1	127mm	127mm	40mm	130
UB 160x6/1	165mm	165mm	6mm	120
UB 160x12/1	164mm	164mm	12mm	90
UB 160x25/1	163mm	163mm	25mm	100
UB 160x40/1	160mm	160mm	40mm	70
UB 200x6/1	208mm	208mm	6mm	80
UB 200x12/1	207mm	207mm	12mm	50
UB 200x25/1	205mm	205mm	25mm	60
UB 200x40/1	200mm	200mm	40mm	40
UB 250x6/1	264mm	264mm	6mm	40
UB 250x12/1	263mm	263mm	12mm	28
UB 250x25/1	260mm	260mm	25mm	35


UniBox rectangular

Article	Length	Width	Height	PU
UB 52x36x5/1	52mm	36mm	5mm	1200
UB 18x55x7/1	55mm	18mm	7mm	1000
UB 18x55x25/1	55mm	18mm	25mm	500
UB 18x69x7/1	69mm	18mm	7mm	1000
UB 18x69x25/1	69mm	18mm	25mm	500
UB 69x45x26/1	69mm	45mm	26mm	350
UB 69x45x41/1	69mm	45mm	41mm	225
UB 31x82x11/1	82mm	31mm	11mm	650
UB 31x82x40/1	82mm	31mm	40mm	330
UB 53x107x18/1	107mm	53mm	18mm	250
UB 53x107x40/1	107mm	53mm	40mm	180
UB 59x126x22/1	126mm	59mm	22mm	140
UB 59x126x40/1	126mm	59mm	40mm	120
UB 82x138x27/1	138mm	82mm	27mm	90
UB 82x138x40/1	138mm	82mm	40mm	60



Multi-talent, all-rounder, universal genius. This box can be used for almost everything

RoseCase ErgoLine

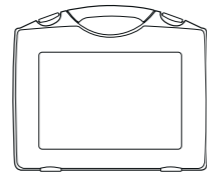

Plastic cases with top functionality and sales boosting design

- modern design gives striking impact
- extremely adaptable to specific market or product requirements through numerous color combinations and various options of interior design – including plain, convoluted, or die-cut foam inserts
- lids and bases can be custom printed or embossed with company logo
- with polished surface label area
- optional transparent lid provides perfect product visibility
- ergonomic handle, optional with soft-touch
- concealed hinges, sturdy, cleverly integrated feet
- excellent for stacking and storage
- numerous secondary uses



 Modern design combined with top functionality – a packaging solution in a class of its own

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.



PACKAGING

Packed assembled



MATERIAL

- Polypropylene (PP)
- Standard version with polished label area



COLORS


- Standard colors "M" & "L": black, red, blue, silver-gray and transparent
- Standard colors "S": transparent
- Closures in black, blue, red, silver-gray or in transparent
- Additional colors and color combinations upon request




ACCESSORIES

- Please inquire for customized foam inserts

INSIDE DEPTHS

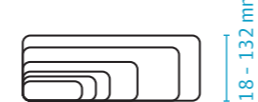

INNER DIMENSIONS

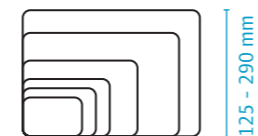

180 - 421 mm

PRODUCTION ENVIRONMENT

- Industrial standard

STANDARD COLORS



Product Range

RoseCase ErgoLine


Article	Outer width	- height	- depth	Inner width	- height	- depth	PU
RCEL 125 S	197mm	167mm	23mm	180mm	125mm	18mm	45
RCEL 125 SM	197mm	167mm	35mm	180mm	125mm	30mm	30
RCEL 125 M	197mm	167mm	47mm	180mm	125mm	42mm	22
RCEL 155 S	244mm	207mm	29mm	225mm	155mm	23mm	30
RCEL 155 SM	244mm	207mm	44mm	225mm	155mm	38mm	20
RCEL 155 M	244mm	207mm	59mm	225mm	155mm	53mm	15
RCEL 175 S	275mm	231mm	31mm	255mm	175mm	26mm	30
RCEL 175 SM	275mm	231mm	48mm	255mm	175mm	43mm	20
RCEL 175 M	275mm	231mm	65mm	255mm	175mm	60mm	15
RCEL 195 S	307mm	260mm	35mm	284mm	195mm	30mm	40
RCEL 195 SM	307mm	260mm	55mm	284mm	195mm	48mm	25
RCEL 195 M	307mm	260mm	74mm	284mm	195mm	66mm	18
RCEL 220 S	346mm	292mm	39mm	320mm	220mm	34mm	28
RCEL 220 SM	346mm	292mm	62mm	320mm	220mm	56mm	18
RCEL 220 M	346mm	292mm	84mm	320mm	220mm	78mm	13
RCEL 250 M	393mm	331mm	95mm	365mm	250mm	90mm	8
RCEL 290 M	456mm	384mm	110mm	421mm	290mm	100mm	5
RCEL 290 ML	456mm	384mm	127mm	421mm	290mm	116mm	4
RCEL 290 L	456mm	384mm	144mm	421mm	290mm	132mm	4

first letter stands for cover, second letter stands for base; e.g. SM = cover size S and base size M

Foam inserts for RoseCase ErgoLine

Article	Suitable for	Height	PU	Article	Suitable for	Height	PU
SE Burl RCEL 125/10	RCEL 125	10mm	1	SE flat foam RCEL 125/10	RCEL 125	10mm	1
SE Burl RCEL 125/15	RCEL 125	15mm	1	SE flat foam RCEL 125/15	RCEL 125	15mm	1
SE Burl RCEL 125/20	RCEL 125	20mm	1	SE flat foam RCEL 125/20	RCEL 125	20mm	1
SE Burl RCEL 155/10	RCEL 155	10mm	1	SE flat foam RCEL 155/10	RCEL 155	10mm	1
SE Burl RCEL 155/15	RCEL 155	15mm	1	SE flat foam RCEL 155/15	RCEL 155	15mm	1
SE Burl RCEL 155/20	RCEL 155	20mm	1	SE flat foam RCEL 155/20	RCEL 155	20mm	1
SE Burl RCEL 155/25	RCEL 155	25mm	1	SE flat foam RCEL 155/25	RCEL 155	25mm	1
SE Burl RCEL 175/10	RCEL 175	10mm	1	SE flat foam RCEL 175/10	RCEL 175	10mm	1
SE Burl RCEL 175/15	RCEL 175	15mm	1	SE flat foam RCEL 175/15	RCEL 175	15mm	1
SE Burl RCEL 175/25	RCEL 175	25mm	1	SE flat foam RCEL 175/30	RCEL 175	30mm	1
SE Burl RCEL 175/30	RCEL 175	30mm	1	SE flat foam RCEL 195/15	RCEL 195	15mm	1
SE Burl RCEL 195/15	RCEL 195	15mm	1	SE flat foam RCEL 195/20	RCEL 195	20mm	1
SE Burl RCEL 195/30	RCEL 195	30mm	1	SE flat foam RCEL 195/30	RCEL 195	30mm	1
SE Burl RCEL 195/35	RCEL 195	35mm	1	SE flat foam RCEL 195/35	RCEL 195	35mm	1
SE Burl RCEL 220/15	RCEL 220	15mm	1	SE flat foam RCEL 220/15	RCEL 220	15mm	1
SE Burl RCEL 220/20	RCEL 220	20mm	1	SE flat foam RCEL 220/25	RCEL 220	25mm	1
SE Burl RCEL 220/30	RCEL 220	30mm	1	SE flat foam RCEL 220/35	RCEL 220	35mm	1
SE Burl RCEL 220/35	RCEL 220	35mm	1	SE flat foam RCEL 220/40	RCEL 220	40mm	1
SE Burl RCEL 220/40	RCEL 220	40mm	1	SE flat foam RCEL 250/25	RCEL 250	25mm	1
SE Burl RCEL 250/40	RCEL 250	40mm	1	SE flat foam RCEL 250/45	RCEL 250	45mm	1
SE Burl RCEL 250/50	RCEL 250	50mm	1	SE flat foam RCEL 290/50	RCEL 290	50mm	1
SE Burl RCEL 290/40	RCEL 290	40mm	1	SE flat foam RCEL 290/65	RCEL 290	65mm	1
SE Burl RCEL 290/50	RCEL 290	50mm	1	SE flat foam RCEL 290/25	RCEL 290 L	25mm	1
SE Burl RCEL 290/60	RCEL 290	60mm	1	SE flat foam RCEL 290/25	RCEL 290 M	25mm	1

RoseCaseProTec

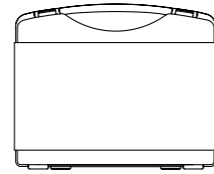

Exclusive plastic case designed to perfectly protect your product

- top quality and technically designed case prevents case shells from bending or twisting
- extremely adaptable to specific market or product requirements through numerous color combinations and various options of interior design – including plain, convoluted, or die-cut foam inserts
- lids and bases can be custom printed or embossed with company logo
- optional transparent lid provides ideal product visibility
- strong and durable 180° opening closures
- concealed hinges and sturdy, cleverly integrated feet
- special location points provide safe stacking and storage
- numerous secondary uses



 **All in one (case) – maximum protection with an appealing product presentation**

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.



PACKAGING

Packed assembled



MATERIAL

- Polypropylene (PP)
- Standard version with polished label area



COLORS

- Standard colors: cover/base in black and transparent; closures in black and red
- RCPT 105: cover with integrated closures; color of closures corresponds to the color of the cover
- Additional colors and color combinations upon request



ACCESSORIES


- Please inquire for customized foam inserts

INNER DEPTHS



35 - 156 mm

INNER DIMENSIONS

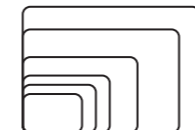

175 - 434 mm

PRODUCTION ENVIRONMENT

- Industrial standard

STANDARD COLORS



Product Range

RoseCase ProTec

Article	Outer width	- height	- depth	Inner width	- height	- depth	PU
RCPT 105 M	187mm	143mm	42mm	175mm	105mm	35mm	28
RCPT 136 M	240mm	183mm	48mm	228mm	136mm	42mm	24
RCPT 161 L	262mm	210mm	72mm	244mm	161mm	66mm	10
RCPT 175 M	275mm	225mm	84mm	261mm	175mm	77mm	9
RCPT 222 S	338mm	273mm	84mm	326mm	222mm	77mm	12
RCPT 222 SM	338mm	273mm	123mm	326mm	222mm	116mm	8
RCPT 222 M	338mm	273mm	163mm	326mm	222mm	156mm	6
RCPT 237 M	395mm	295mm	107mm	382mm	237mm	100mm	5
RCPT 296 M	456mm	357mm	140mm	434mm	296mm	132mm	4

Foam inserts for RoseCase ProTec

Article	Shape	Suitable for	Outer height	PU
SE Burl RCPT 105/10	convoluted foam	RCPT 105	10mm	1
SE Burl RCPT 105/15	convoluted foam	RCPT 105	15mm	1
SE Burl RCPT 105/20	convoluted foam	RCPT 105	20mm	1
SE Burl RCPT 136/20	convoluted foam	RCPT 136	20mm	1
SE Burl RCPT 161/30	convoluted foam	RCPT 161	30mm	1
SE Burl RCPT 175/35	convoluted foam	RCPT 175	35mm	1
SE Burl RCPT 175/40	convoluted foam	RCPT 175	40mm	1
SE Burl RCPT 222/25	convoluted foam	RCPT 222	25mm	1
SE Burl RCPT 222/35	convoluted foam	RCPT 222	35mm	1
SE Burl RCPT 222/40	convoluted foam	RCPT 222	40mm	1
SE Burl RCPT 222/50	convoluted foam	RCPT 222	50mm	1
SE Burl RCPT 222/60	convoluted foam	RCPT 222	50mm	1
SE Burl RCPT 237/50	convoluted foam	RCPT 237	50mm	1
SE Burl RCPT 296/25	convoluted foam	RCPT 296	25mm	1
SE Burl RCPT 296/40	convoluted foam	RCPT 296	40mm	1
SE Burl RCPT 296/60	convoluted foam	RCPT 296	60mm	1
SE Flat foam RCPT 105/10	flat foam	RCPT 105	10mm	1
SE Flat foam RCPT 105/20	flat foam	RCPT 105	20mm	1
SE Flat foam RCPT 136/15	flat foam	RCPT 136	15mm	1
SE Flat foam RCPT 136/20	flat foam	RCPT 136	20mm	1
SE Flat foam RCPT 161/15	flat foam	RCPT 161	15mm	1
SE Flat foam RCPT 161/30	flat foam	RCPT 161	30mm	1
SE Flat foam RCPT 175/25	flat foam	RCPT 175	25mm	1
SE Flat foam RCPT 175/35	flat foam	RCPT 175	35mm	1
SE Flat foam RCPT 222/20	flat foam	RCPT 222	20mm	1
SE Flat foam RCPT 222/35	flat foam	RCPT 222	35mm	1
SE Flat foam RCPT 222/50	flat foam	RCPT 222	50mm	1
SE Flat foam RCPT 222/75	flat foam	RCPT 222	75mm	1
SE Flat foam RCPT 237/35	flat foam	RCPT 237	35mm	1
SE Flat foam RCPT 237/50	flat foam	RCPT 237	50mm	1
SE Flat foam RCPT 296/25	flat foam	RCPT 296	25mm	1
SE Flat foam RCPT 296/40	flat foam	RCPT 296	40mm	1
SE Flat foam RCPT 296/60	flat foam	RCPT 296	60mm	1

AQRADO®

Versatile drawer storage system – modular, customized and installed in seconds!

- wide product range to suit individual applications
- suitable for every drawer size
- neat arrangement means quick access and efficient working
- perfect protection of instruments and tools
- numerous modules for universal applications as well as for special contents such as laboratory equipment, stationery...
- base pad can be cut in all sizes in 24 mm increment spacings
- solvent resistant



PLEASE INQUIRE FOR INDIVIDUAL DRAWER SETS IN SMALL QUANTITIES



PACKAGING

Parts are packed separately in individual cartons



MATERIAL


- Impact resistant ABS plastics




COLORS

- Standard color black and gray

OUTER DEPTHS

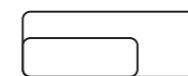


OUTER DIMENSIONS

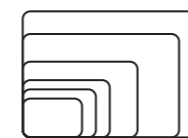

PRODUCTION ENVIRONMENT

- Industrial standard

STANDARD COLORS

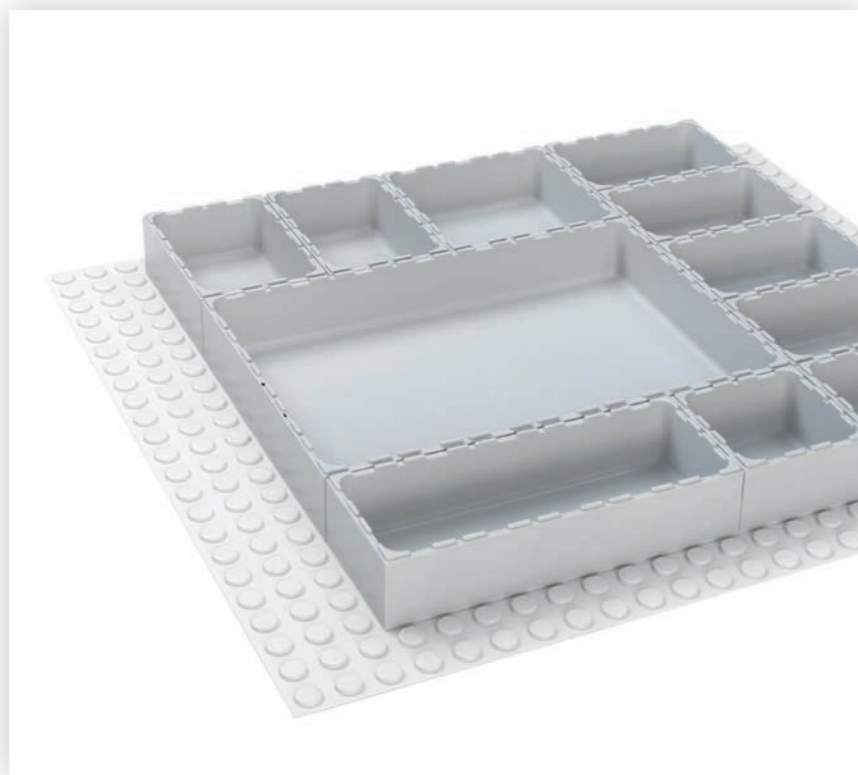

Compartmentalized thinking can sometimes make sense - if it ensures perfect organization

The proof of the suitability of the packaging in combination with the medical device to be packed must be provided by the manufacturer.

Product Range


Aqrado® Sets

Article	Description	Outer length	Outer width	Outer height	PU
AQ-Set 5443	12 modules + base pad	600mm	408mm	24mm	1
AQ-Set 5444	11 modules + base pad	600mm	576mm	48mm	1
AQ-Set 5445	20 modules + base pad	600mm	600mm	48mm	1
AQ-Set 5446	24 modules + base pad	528mm	480mm	24mm	1

Aqrado® Modules

Article	Description	Outer length	width	height	PU
AQ-Multi 0102	3 compartments	192mm	96mm	24mm	84
AQ-Multi 0103	2 compartments	288mm	96mm	24mm	56
AQ-Multi 0113	8 compartments	192mm	144mm	24mm	56
AQ-Uni 0100		96mm	48mm	24mm	120
AQ-Uni 0101		96mm	96mm	24mm	60
AQ-Uni 0104		192mm	96mm	24mm	84
AQ-Uni 0105		96mm	96mm	48mm	84
AQ-Uni 0106		192mm	96mm	48mm	72
AQ-Uni 0107		240mm	96mm	48mm	57
AQ-Uni 0108		288mm	96mm	48mm	48
AQ-Uni 0109		336mm	96mm	48mm	39
AQ-Uni 0110		288mm	96mm	24mm	56
AQ-Uni 0111		144mm	96mm	48mm	56
AQ-Uni 0112		144mm	96mm	24mm	40
AQ-Uni 0114		384mm	144mm	48mm	24
AQ-Uni 0115		336mm	240mm	48mm	15
AQ-Uni 0116		144mm	144mm	48mm	64
AQ-Uni 0117		192mm	144mm	48mm	48
AQ-Uni 0118	module for 5 dividers	144mm	144mm	48mm	64
AQ-Uni 0119	module for 7 dividers	192mm	144mm	48mm	48
AQ-Uni 0120	module for 15 dividers	384mm	144mm	48mm	24
AQ-Uni 0121		192mm	72mm	48mm	56
AQ-Uni 0122		192mm	72mm	24mm	40
AQ-Vario 1114	end-piece, 8 compartments	192mm	48mm	24mm	60
AQ-Vario 1115	end-piece, 8 compartments	192mm	73mm	24mm	40
AQ-Vario 1116	end-piece, 3 compartments	96mm	96mm	24mm	60
AQ-Vario 1117	end-piece, 2 compartments	96mm	96mm	24mm	60
AQ-Vario 1118	end-piece, 1 compartments	73mm	48mm	24mm	80
AQ-Vario 1119	end-piece, 1 compartments	96mm	96mm	24mm	60
AQ-Vario 1121	end-piece, 1 compartments	92mm	72mm	48mm	32
AQ-Vario 1122	end-piece, 1 compartments	92mm	72mm	24mm	80

Vario modules: two parts needed

Accessories see next page 



Product Range

Aqurado® accessories

Article	Description	Outer length	Outer width	Outer height	PU
AQ-Extra 4100	cutting guide/aid				25
AQ-Extra 4101	multi-purpose knife				60
AQ-Extra 4102	adhesive pads (6 pcs.)				25
AQ-Extra 4103	adhesive pads (34 pcs.)				50
AQ-Extra 4104	label tab (vertical)	45mm	13mm		50
AQ-Extra 4105	label tab (horizontal)	48mm	15mm		50
AQ-Extra 4107	eccentric fixing pins				20
AQ-Extra 4390	assembly kit				10
AQ-Extra 4400	cover	96mm	96mm		10
AQ-Extra 4401	cover	96mm	48mm		10
AQ-Extra 4402	cover	144mm	96mm		10
AQ-Extra 4403	cover	192mm	96mm		10
AQ-Extra 4500	divider			60mm	10
AQ-Extra 4501	divider			100mm	10

Aqurado® base pad

Article	Outer length	Outer width	PU
AQ-Base 3907	1584mm	624mm	16



PACKAGING

Parts are packed separately in individual cartons



MATERIAL

- Impact resistant ABS plastics



COLORS


- Standard color
- black and gray

INNER DEPTHS



24 and 48 mm

INNER DIMENSIONS

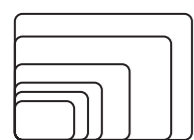

48 - 240 mm

73 - 384 mm

PRODUCTION ENVIRONMENT

- Industrial standard

STANDARD COLORS





HEADQUARTERS
Germany

Production locations

GERMANY

rose plastic medical packaging GmbH

Rupolzer Straße 30
88138 Hergensweiler/Lindau
Tel. +49 8388 92345-0
info@rose-medipack.de

USA

rose plastic medical packaging USA LLLP

P.O. Box 210, California, PA 15419
545 Technology Drive, Coal Center, PA 15423
Tel. +1 724 938 8530
info@rose-medipack.us

CHINA

rose plastic (Kunshan) Co., Ltd.

Donghui Road No. 101, Zhoushi Town,
Kunshan, Jiangsu Province
Tel. +86 512 5766 7700
info@rose-plastic.cn

BRAZIL

rose plastic Brasil Embalagens Plásticas Ltda.

Av. Garabed Ganianiam 514, Jardim Topázio
CEP 18087-340 Sorocaba SP
Tel. +55 15 3142 3050
info@rose-plastic.com.br

INDIA

rose plastic India Pvt. Ltd.

Pargaon, Khandala, Satara - 412802
Maharashtra - India
Tel. +91 21 6929 9151
info@rose-plastic.in

Subsidiaries

UNITED KINGDOM

rose plastic UK Ltd.

Unit 4, Bessemer Way, Bessemer
Business Park, Rotherham, S60 1EN
Tel. +44 1709 7217 94
info@rose-plastic.co.uk

FRANCE

rose plastic France S.A.R.L.

Parc d'Activités du Rotey
73460 Notre Dame des Millières
Tel. +33 479 3848 01
info@rose-plastic.fr

SPAIN

rose plastic Iberia, S.L.U.

P.A.E. Ibarra Barri, Pab. 15B
48940 Leioa (Vizcaya)
Tel. +34 94 480 66 61
info@rose-plastic.es

REPUBLIC OF KOREA

rose plastic Korea Co., Ltd.

Bucheon Techno Park 402-B01, B02,
Pyung-Cheon Ro 655, Won-mi-gu,
Bucheon City, Gyung-gi-do, ROK 14502
Tel. +82 32 324-8332
info@rose-plastic.kr

ITALY

rose plastic Italia S.r.l.

Viale De Gasperi 103
20017 Rho (MI)
Tel. +39 02 9390 9450
info@rose-plastic.it



10/2024